

Surgical technique


balanSys BICONDYLAR

LIS 4in1

SMarT instruments

For healthcare professional use only. The illustrated image does not represent a connection between the use of the medical device described, nor its performance.

Preservation in motion

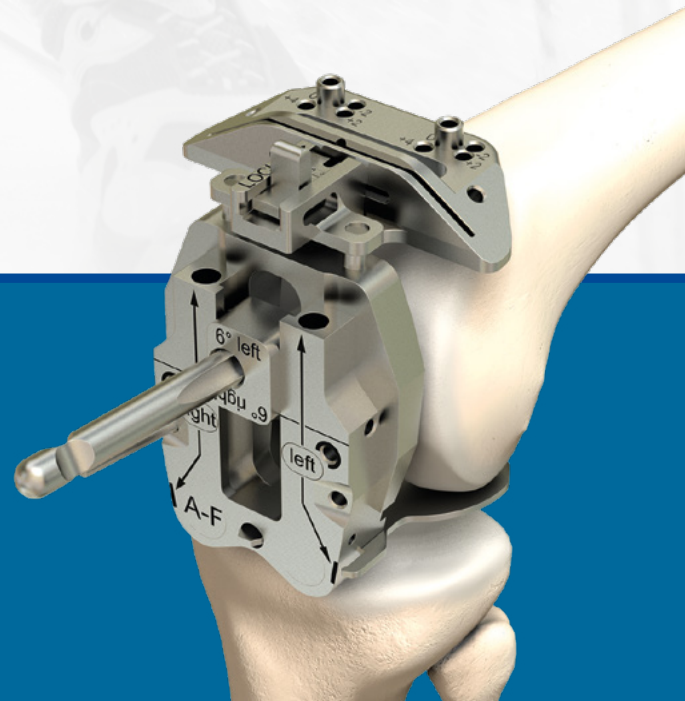

*Building on our heritage
Moving technology forward
Step by step with our clinical partners
Towards a goal of preserving mobility*



Preservation in motion

As a Swiss company, Mathys is committed to this guiding principle and pursues a product portfolio with the goal of further developing traditional philosophies with respect to materials or design in order to address existing clinical challenges. This is reflected in our imagery: traditional Swiss activities in conjunction with continuously evolving sporting equipment.

Table of contents

Introduction	4
Indications and contraindications	5
Preoperative planning	6
Aim of the intervention and surgical approach	7
Patient preparation	7
Overview of the surgical technique	8
Surgical technique	10
1. Femoral osteotomy	10
2. Tibial osteotomy	18
3. Preparation and implantation	22
3.1 balanSys CR, UC and RP	22
3.2 balanSys PS	27
Appendix	33
1 – Intramedullary tibial osteotomy	33
2 – balanSys 3-Peg Patella	35
3 – Item numbers of the balanSys implants	38
4 – balanSys implant size and compatibility	50
5 – Item numbers of the balanSys instruments	51
6 – Item numbers of the balanSys measuring templates	80
Symbols and abbreviations	81

Remark

Please make yourself familiar with the handling of the instruments, the product-related surgical technique and the warnings, the safety notes as well as the recommendations of the instruction leaflet before using an implant manufactured by Mathys Ltd Bettlach. Make use of the Mathys user training and proceed according to the recommended surgical technique.

Introduction

The balanSys BICONDYLAR system is a modular joint replacement system for knees.

A balanSys BICONDYLAR system consists of a cemented or non-cemented femoral component, a cemented symmetrical tibial tray and a tibial inlay. A cemented patellar component is optional. For femoral and tibial metal components a TiNbN coated option is available. The tibial inlay range includes one mobile bearing configuration (Rotating Platform (RP)) and three different fixed bearing configurations (Cruciate Retaining (CR), Ultra Congruent (UC) and Posterior Stabilized (PS)).

For tibial inlays two materials can be chosen: Standard UHMWPE or vitamys, the vitamin E stabilized PE, which possess a high wear resistance and is well protected against oxidation.¹

Since 1997, the balanSys BICONDYLAR system has proven its clinical worth with good clinical results.²

With the balanSys knee system, Mathys Ltd Bettlach offers a wide range of components that correspond with the patient's anatomical conditions and functional requirements of the knee joint.

The balanSys implants, instruments and surgical techniques are designed to meet the surgeons' demands on prostheses with respect to kinematics, ligament balancing, stability and long-term survival.²

The high quality instruments and implants are made for continuously excellent clinical results. Quality made in Switzerland.

¹ Data on file Mathys Ltd Bettlach.

² Superior long-term survival for fixed bearing compared with mobile bearing in ligament-balanced total knee arthroplasty. Heesterbeek, P.J.C., van Houten, A.H., Klenk, J.S. et al. Knee Surg Sports Traumatol Arthrosc, online 07 April 2017.

Indications and contraindications

Indications

- Painful and / or disabling joint disease of the knee resulting from osteoarthritis, avascular necrosis, inflammatory arthritis or post-traumatic arthritis
- Revision of previous knee replacement

Contraindications

- Local or general infection
- Any soft tissue, ligament, nerve or vessel insufficiency which may create an unacceptable risk of prosthesis instability, prosthesis fixation failure and/or complications in post-operative care
- Compromised bone stock due to bone loss or bone defects and/or insufficient bone substance, which cannot provide adequate support and/or fixation for the prosthesis
- Hypersensitivity to materials used
- Skeletal immaturity
- Genu recurvatum
- Insufficiency of the extensor mechanism
- Progressive neoplastic disease

For further information, please refer to the instructions for use or ask your Mathys representative.

Preoperative planning

Preoperative planning aims

- Identification and quantification of deformities and bone defects as well as the osteophytes
- Determination of the angle between the anatomical axis and the mechanical axis
- Preliminary size determination of the femoral and tibial prosthesis using planning templates

The preoperative planning includes

- Anteroposterior (AP) and lateral X-rays of the knee joint and tangential X-rays of the patella
- Full leg series
- Total leg planning template
- X-ray template

Aim of the intervention and surgical approach

- Intraoperative correction of axial deviations in the frontal plane of the leg along the mechanical axis, whereby the joint line should be at a right angle with this axis
- Reconstruction of the physiological axis ratios
- Prosthesis-related kinematics
 - physiological joint line
 - sufficient medial and lateral stability in extension and flexion
 - correctly centred and balanced patellofemoral joint
 - free movement: from maximal extension to maximal possible flexion

The choice of approach is dependent on the axial malposition (varus/valgus).

Patient preparation

The surgery is carried out on patients under general or spinal anaesthesia, while an adequate muscle relaxation is being required.

Postoperative pain is reduced without the use of a tourniquet. If it is necessary to apply a tourniquet, it should be placed on the proximal thigh and inflated with the knee in hyperflexion. That will keep the maximal portion of the quadriceps below the level of the tourniquet.


Place the patient in supine position.

Flex the knee into a 90° position.




Use a supporting roll on the table and a lateral support to facilitate extension and flexion of the leg.

Overview of the surgical technique

1. Femoral osteotomy

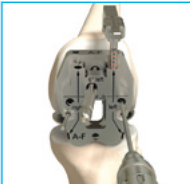				Opening of the intramedullary canal and insertion of the Intramedullary Rod. Attachment of the Femoral Sizing/Rotation Guide LIS. Determination of rotational orientation and femoral size. > Page 10
				Drilling of two holes for the 4in1 Cutting Block. Fixation of the Distal Femoral Cutting Block. Distal Osteotomy. > Page 14
			Insertion of the 4in1 Cutting Block. Control of the planned resection depth. Anterior and posterior femoral osteotomies with the chamfer cuts. > Page 16	

2. Tibial osteotomy

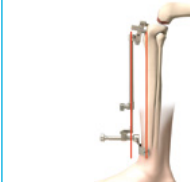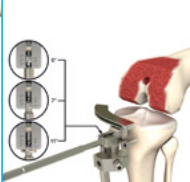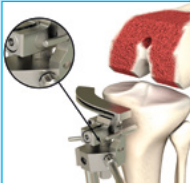				Application of the Tibial Reference System parallel to the anterior tibial cortex and alignment. Adjustment of the posterior slope. Determination of the joint line and fixation of the Tibial Reference System. > Page 18
				Adjustment of the resection depth. Tibial osteotomy. Determination of the tibial plateau size. Remark Place retractors to protect ligaments during tibial resection. > Page 19

3.1 Preparation and implantation: balanSys CR, UC and RP

				Preparation of trochlea groove. Insertion of Tibial Template and Trial PE Inlay. Insertion of Trial Femur. Trial reduction of the knee joint. > Page 22
				I. Preparation of the femoral anchor pins II. Preparation of the tibial medullary space III. Preparation of the fin > Page 24
				I. Insertion of balanSys Tibial Plateau II. Insertion of the balanSys Inlay III. Insertion of the balanSys Femur > Page 25

3.2 Preparation and implantation: balanSys PS

				Preparation of femoral box. Insertion of the Trial Femur. Insertion of the Tibial Template and PE Trial Inlay. Trial reduction of the knee joint. > Page 28
				I. Preparation of the tibial medullary space II. Preparation of the fin > Page 31
				I. Insertion of the balanSys PS Tibial Plateau II. Insertion of the balanSys PS Femur III. Insertion of the balanSys PS Inlay > Page 32

Remark

1.27mm sawblades are compatible with balanSys instruments. For sawblades distributed by Mathys refer to the brochure 336.030.032 «Sterile Sawblades».

Surgical technique

1. Femoral osteotomy

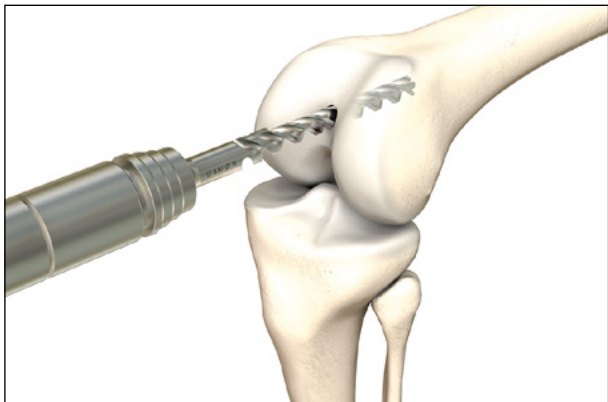


Fig. 1 Opening of the medullary cavity

Intramedullary alignment of the femoral resection

Open the medullary cavity with the 8.5mm Drill Bit (Fig. 1).

Remark

The femoral medullary cavity must be opened in the correct axial extension of its course to avoid incorrect implantation in flexion and extension.

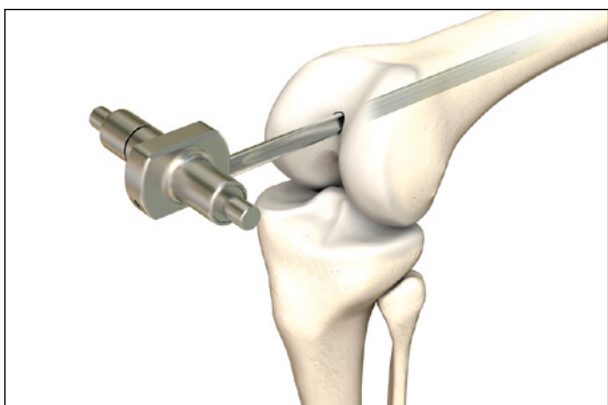


Fig. 2 Introduction of the Intramedullary Rod

Introduce the Intramedullary Rod using the Tibia Reference System Handle. Existing osteophytes must be removed first (Fig. 2).

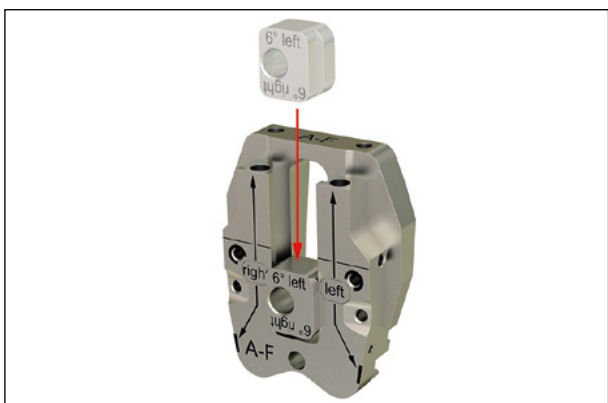


Fig. 3 Introduction of the Axis Alignment Block LIS

Positioning of the Femoral Sizing / Rotation Guide

Select size of Femoral Sizing/Rotation Guide (XS–D or A–F). Introduction of the Axis Alignment Block LIS (3°, 5°, 6°, 7° or 9°) into the Femoral Sizing/Rotation guide (Fig. 3). The angle of the Axis Alignment Block corresponds to the pre-operatively determined angle between the mechanical and anatomical femoral axis (femoral valgus angle).

Remark

Introduce Axis Alignment Block according to right or left leg. «Left» means a left knee; «right» means a right knee. The orientation of the lettering on the Axis Alignment Block must be identical to that on the Femoral Sizing/Rotation Guide.

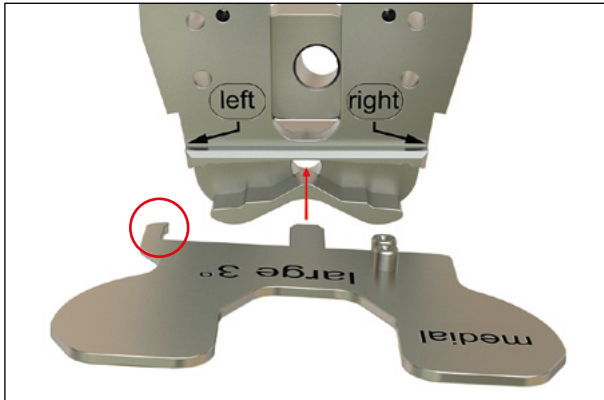


Fig. 4 Introduction of the Femur Rotation Bearing

For rotation adjustment, 3 angles are available. Depending on the angle, the Femur Rotation Bearing 0°, 3° or 5° can be selected.

Introduction of the desired Femur Rotation Bearing small or large (Fig. 4).

Remark

The pointer on the Femur Rotation Bearing must point to the respectively specified side (right or left) on the Femoral Sizing/Rotation guide. For a left knee, the pointer must point to «left», for a right knee, to «right» (Figs. 4 & 5).

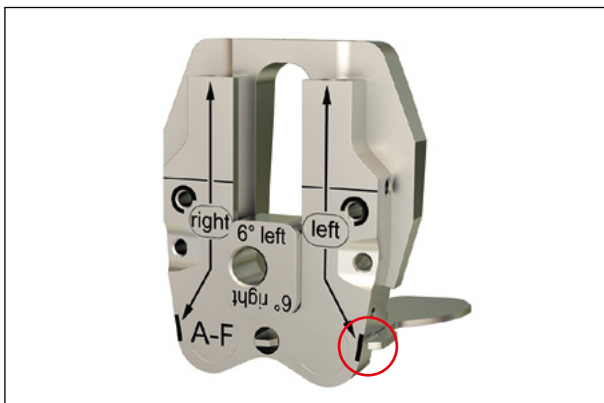


Fig. 5 Position of the pointer

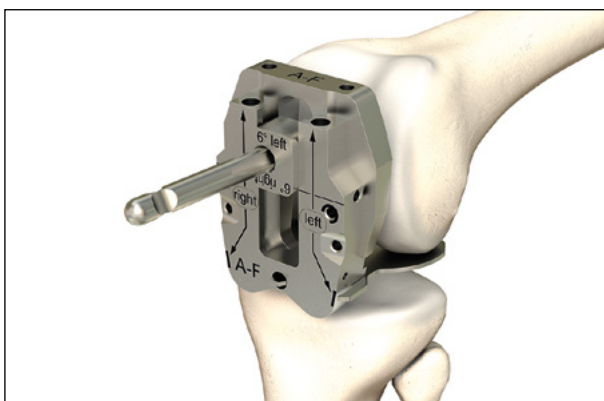


Fig. 6

Push up the Femoral Sizing/Rotation Guide over the Intramedullary Rod through the borehole in the Axis Alignment Block until the instrument stops at a distal condyle (Fig. 6).


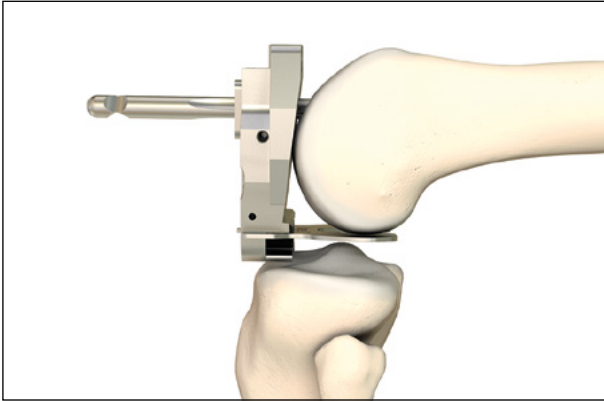


Fig. 7

Push instrument in anterior direction until the Femur Rotation Bearing touches the posterior condyles on both sides (Fig. 7).

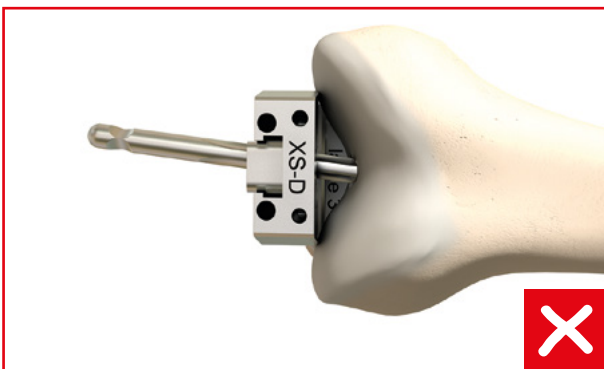

If the Femoral Sizing/Rotation Guide XS–D slides between the condyles in the notch, the Femoral Sizing/Rotation guide A–F must be selected (Fig. 8).

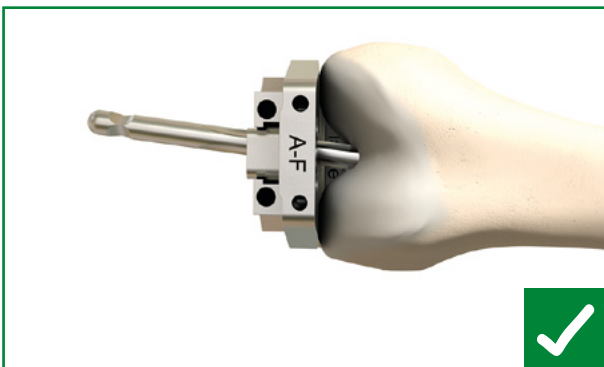


Fig. 8

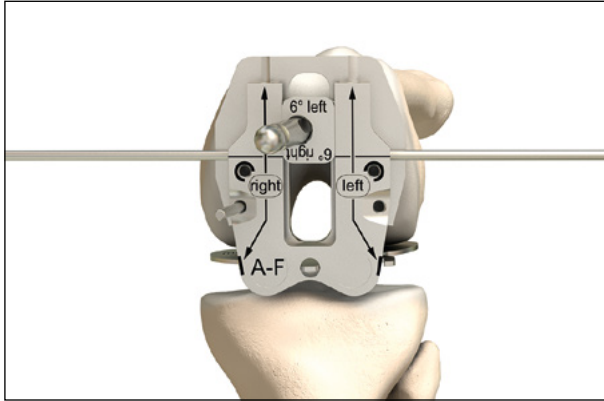


Fig. 9

The Femoral Sizing/Rotation Guide can be accordingly adjusted using the scribed line and the inserted Pins for line extension to the transepicondylar axis (Fig. 9).

After adjustment, the Femoral Sizing/Rotation Guide is fixed with a Pin through the oblique holes in the lower instrument area. The hole is pre-drilled using the 3.2 mm Drill Bit.

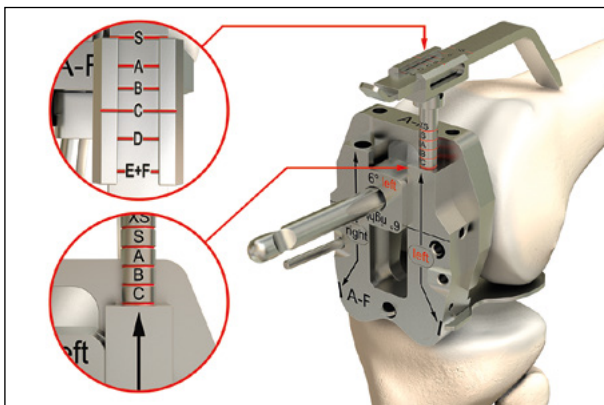


Fig. 10 Determination of the femoral sizes

The femoral size is determined using the Size Feeler LIS, which is placed on the highest point of the femoral metaphysis.

Introduce the Size Feeler LIS on the relevant side (right/ left) of the Femoral Sizing/Rotation Guide (Fig. 10).

The size of the femur is read on the anterior and distal scale of the Size Feeler LIS.

1. Read the marking on the distal scale.
2. Transfer the value read to the anterior scale.
3. The two sizes must match.

Remark

The size of the femur must match the tibial size (Appendix 4).

Remark

The Femoral Feeler LIS proximally touches the position at which the tip of the anterior femoral plate comes to rest after implantation.

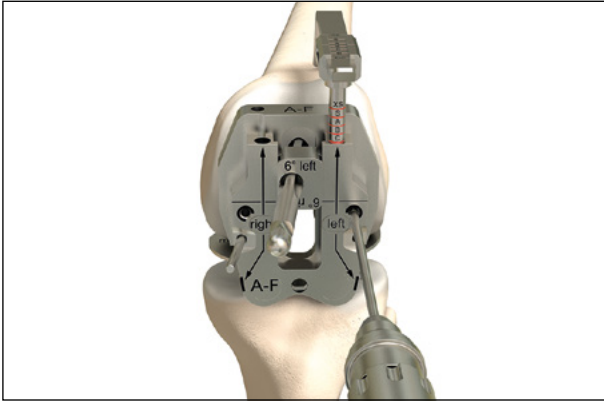


Fig. 11 Drilling of the two receiving holes with the 3.2 mm Drill Bit

After the sizing of the femur, the two holes for receiving the 4in1 Cutting Block are drilled (Fig. 11).

Remark

The two bores of the receiving holes must be deep enough to remain visible even after the distal femoral resection.

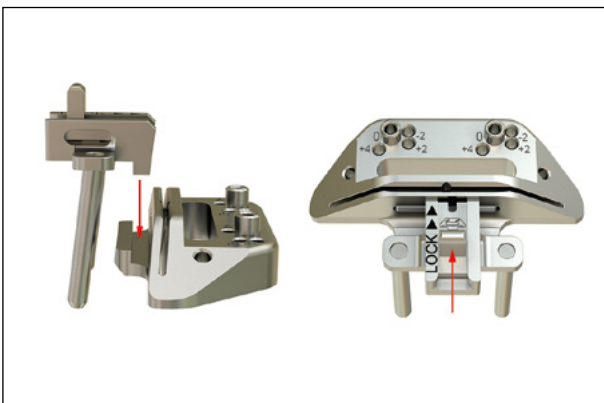


Fig. 12

Place the Distal Cutting Block onto the Holder LIS.

Push the slider in the Holder LIS forward until it latches into place (Fig. 12).

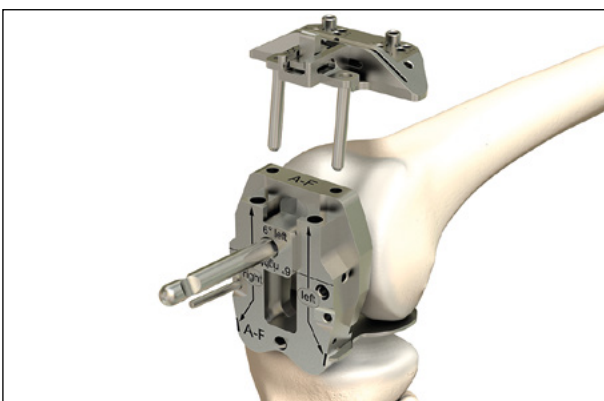


Fig. 13

Introduce the Holder LIS with the Distal Cutting Block into the borehole of the Femoral Sizing/Rotation Guide and push it onto the anterior femur (Fig. 13).

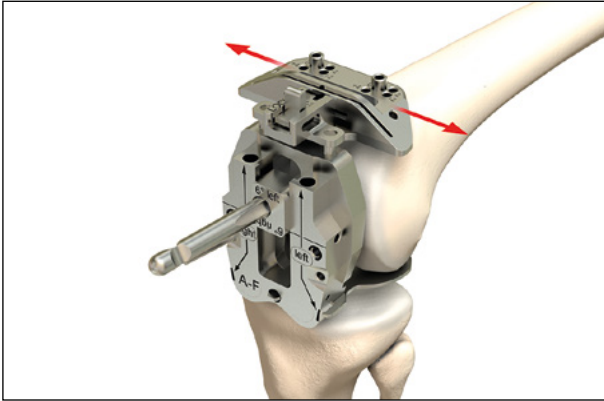


Fig. 14 Positioning the Distal Cutting Block

The Distal Cutting Block can be moved mediolaterally on the basis of the anatomical situation (Fig. 14).

The Distal Cutting Block must butt against the anterior condyle.

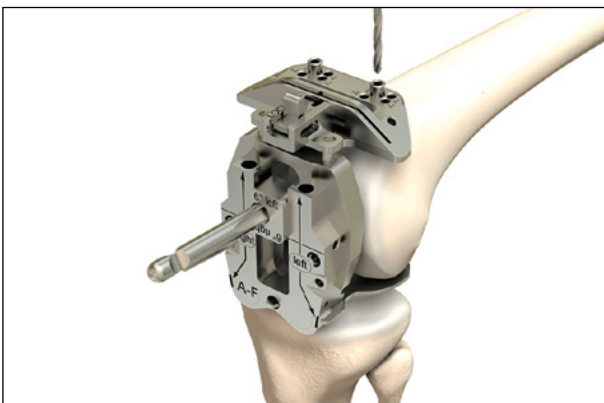


Fig. 15 Pre-drilling and setting of 2 Pins

Pre-drilling of the two reference holes for the Distal Cutting Block. The 0-mark must be paid attention to. Insert two Pins into the anteriorly pre-drilled holes (Fig. 15).

Removal of the instruments in the following order:

1. Loosen the slide on the Holder LIS, then pull out the Holder from the LIS Femoral Sizing/Rotation Guide.
2. Pull out Fixation Pin in the Femoral Sizing/Rotation Guide.
3. Remove Femoral Sizing/Rotation Guide incl. Axis Alignment Block and Femur Rotation Bearing.
4. Pull out Intramedullary Rod.

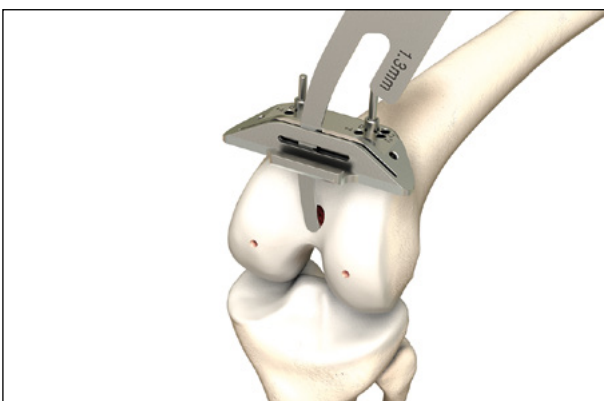


Fig. 16 Control of the planned distal resection surface

Check the planned distal resection plane using the Feeler Plate (Fig. 16).

If necessary, move the Distal Cutting Block by +2/+4/-2 mm.

Remark

The Distal Cutting Block can be moved to obtain constant ligament tension in flexion and extension.

Remark

To ensure optimum stability, the Distal Cutting Block should be fixed with an additional oblique pin.

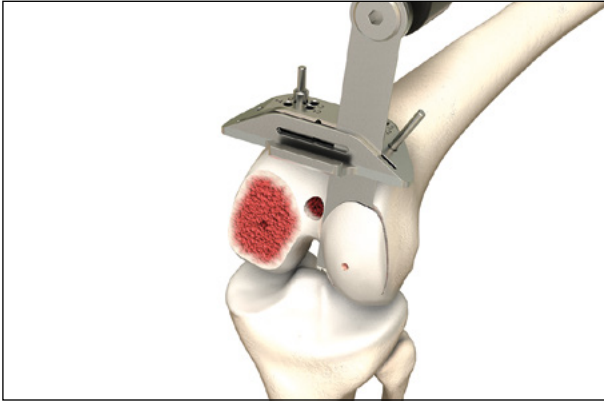


Fig. 17 Performing the distal femoral resection

Performing the distal femoral resection (Fig. 17).

Removal of the Pins and the Distal Cutting Block.

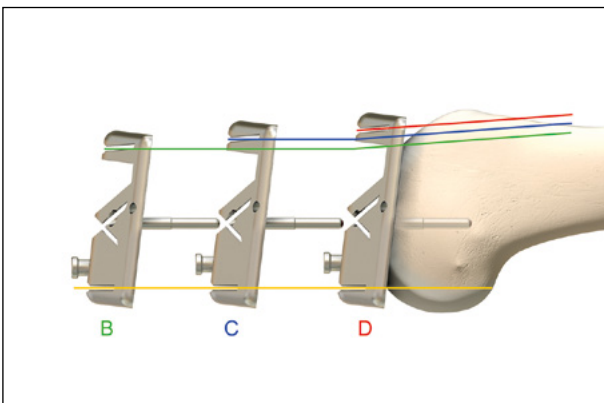


Fig. 18 Insertion of the 4in1 Cutting Block

Remark

The balanSys Femur components **differ 3mm in AP dimension** from one size to the other. The distance between the posterior cut and the pin fixation is exactly the same throughout the entire range of Cutting Blocks, thus leaving the flexion gap constant between sizes (Fig. 18).

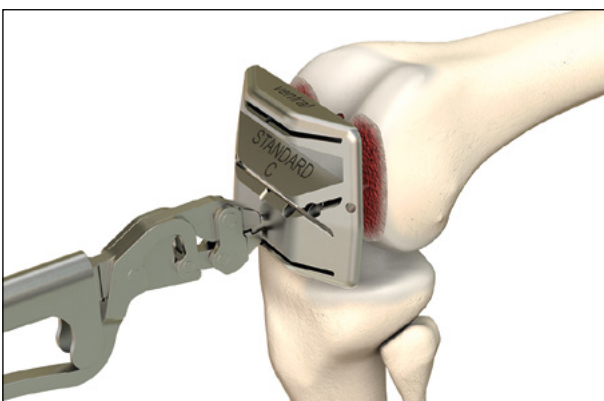


Fig. 19 Insertion of the 4in1 Cutting Block

Execution of the anterior and posterior femoral osteotomies with the chamfer cuts

Place the selected 4in1 Cutting Block in the two pre-drilled holes with the aid of the Pliers until it abuts flat on the distal osteotomy surface (Fig. 19).

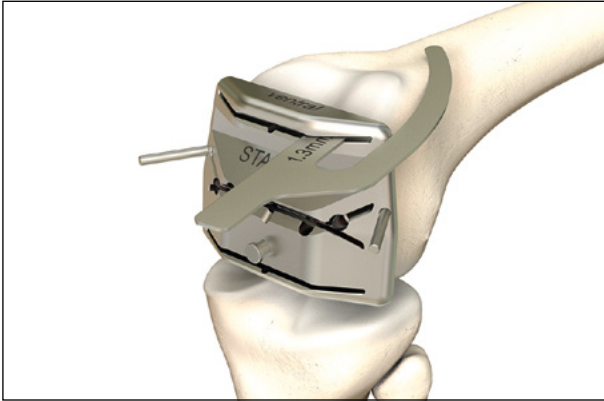


Fig. 20

Check the planned resection plane with the Reference Plate (Fig. 20). Mediolateral fixation of the 4in1 Cutting Block with two Pins for stabilization.

Remark

The 4in1 Cutting Block can be shifted 1.5 mm anterior and 1.5mm posterior. Pre-drill the corresponding holes through the 4in1 Cutting Block. Remove the 4in1 Cutting Block with the aid of the Pliers and place it in the new pre-drilled holes until it abuts flat on the distal osteotomy surface.

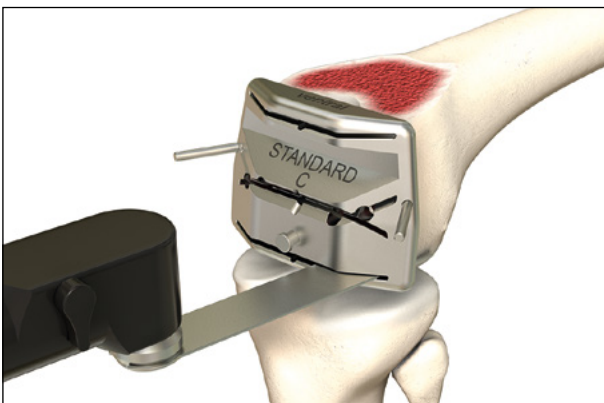


Fig. 21 Execution of the osteotomies

Execution of the osteotomies (Fig. 21) in the following order:

1. Anterior osteotomy
2. Posterior osteotomy
3. Chamfer cuts

Remark

The posterior osteotomies should be performed with the knee at 90° flexion, as this avoids touching the tibial surface with the saw blade and it moves the posterior soft tissues away from the posterior condyles.

Removal of the 4in1 Cutting Block with the aid of the Pliers.

Removal of all marginal bones and osteophytes, particularly in the region of the posterior condyles.

Surgical technique

2. Tibial osteotomy

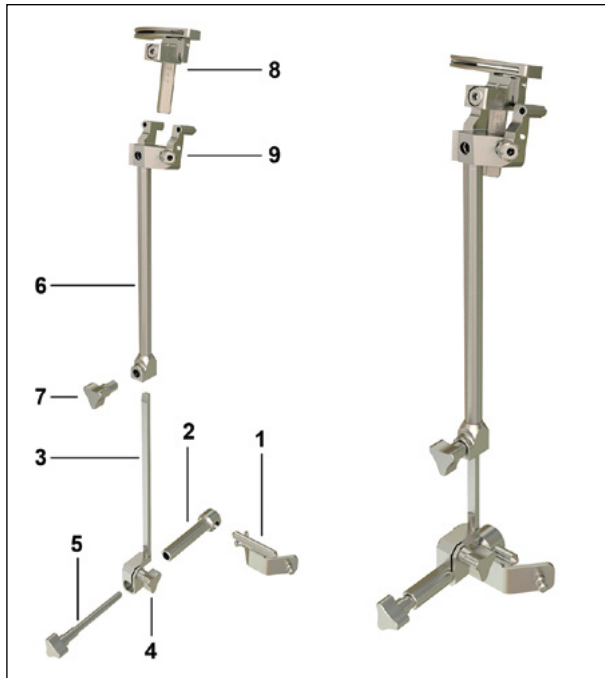


Fig. 22

Assembly of the Extramedullary Tibial Alignment Reference System (Fig. 22)

Move the Ankle Holder (1) onto the Distal Connector (2). Move the Distal Target Device (3) onto the Distal Connector (2) and fix the Ankle Holder with the Tibial Locking Bolt (5). Leave the Connecting Screw (4) slightly open. Join the Distal Target Device (3) with the Proximal Target Device (6) and affix it with the Connecting Screw (7). Insert the Tibial Cutting Guide (8) onto the Proximal Targeting Device (6). Set the Tibial Cutting Guide to zero and tighten the screw (9) with the Hexagonal Screwdriver.

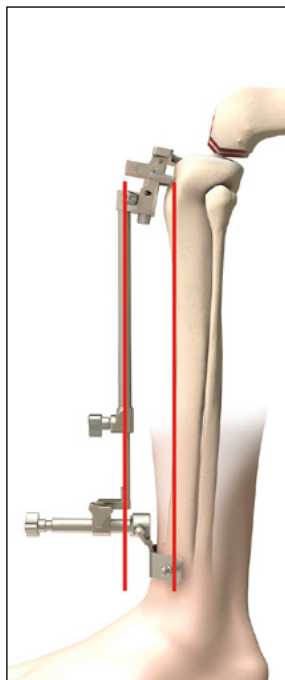


Fig. 23 Application of the Tibial Reference System



Fig. 24 Alignment of the Tibial Reference System

Positioning of the extramedullary Tibial Alignment Reference System

The Tibial Reference System is placed parallel to the anterior tibial cortex (Fig. 23). Align the Tibial Reference System distally with the second metatarsophalangeal bone and fix it with the Silicone Band (Fig. 24). For correct rotation alignment, the centre of the Tibial Cutting Guide must be aligned to the transition of the medial to the middle third of the tibial tuberosity. Alternatively, the Tibial Cutting Guide can be aligned parallel to the line between the posterior edges of the medial and lateral tibial plateau.

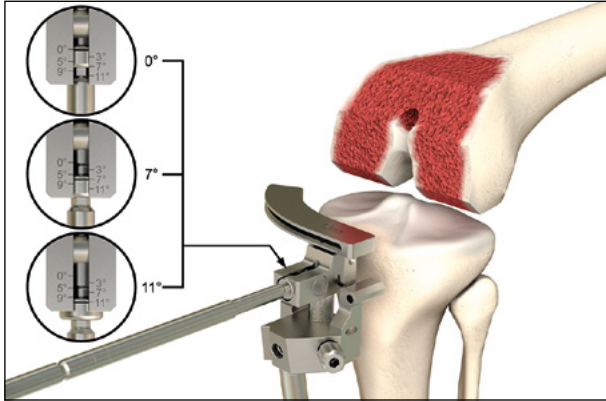


Fig. 25 Adjustment of the dorsal slope

Setting of the posterior slope

Use the Hexagonal Screwdriver to set the posterior slope, according to the anatomy (tactile plate parallel to the best preserved tibial joint surface) (Fig. 25).

Remark

The authors recommend a posterior slope of 7° for a cruciate-retaining implant and up to 5° for a PCL-substituting implant.

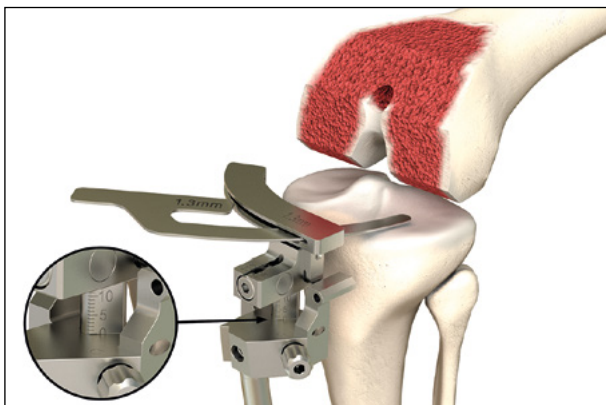


Fig. 26 Setting the original joint line

Determine the original joint line at the level of the best preserved tibial joint surface with the Reference Plate or the Tibial Stylus (Fig. 26).



Fig. 27 Fixation of the Tibial Reference System

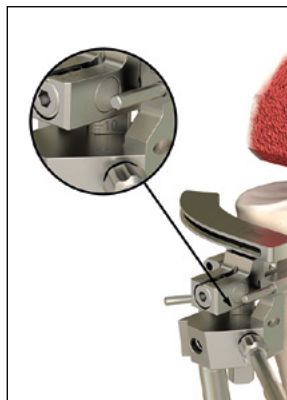


Fig. 28 Adjustment of the resection depth

Fixation of the Tibial Reference System

Affix the Tibial Reference System proximally with at least two Pins in the specified holes (oblique and straight) (Fig. 27). The holes are pre-drilled with the 3.2 mm drill.

Tibial osteotomy

Move the Tibial Cutting Guide 6–8 mm distally with the help of the Hexagonal Screwdriver (Fig. 28).

Check the adjusted osteotomy level with the Reference Plate.

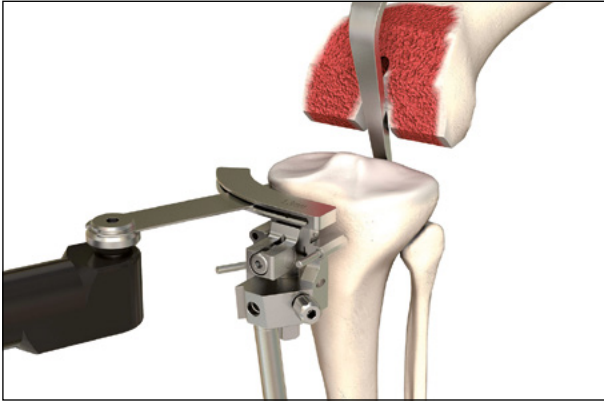


Fig. 29 Execution of the tibial osteotomy

Execution of the tibial osteotomy (Fig. 29) and removal of the posterior meniscus remnant.

Remark

Place retractors to protect ligaments during tibial resection.

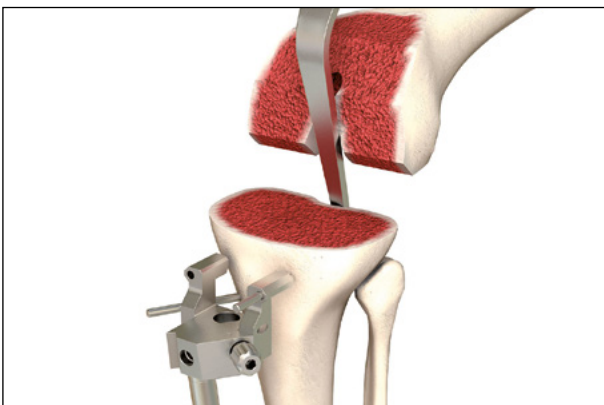


Fig. 30

Remove the Tibial Cutting Guide (Fig. 30).

Remark

If the Tibial Reference System is removed, a Pin should be left in place for future reference in case an additional tibial osteotomy is required.

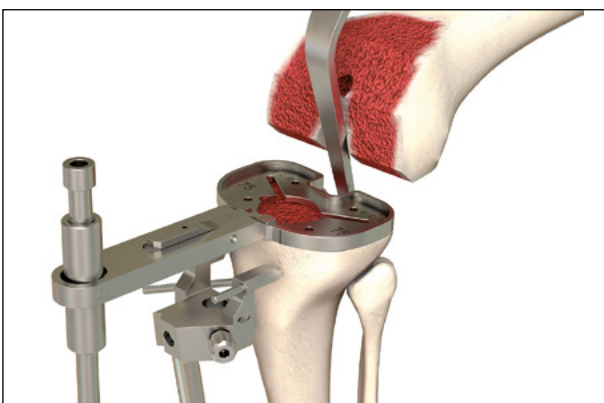


Fig. 31 Determining the tibial prosthesis size

Determining the size of the tibial implant

Determine the tibial prosthesis size with the Tibial Template, taking the rotational alignment into account (Fig. 31). The aim is to provide maximum coverage of the osteotomy surface. Use the Setting Rod to check the axis and the slope of the cutting plane.

In addition to the extramedullary aligned tibial osteotomy, intramedullary alignment is also possible (see Appendix 1).

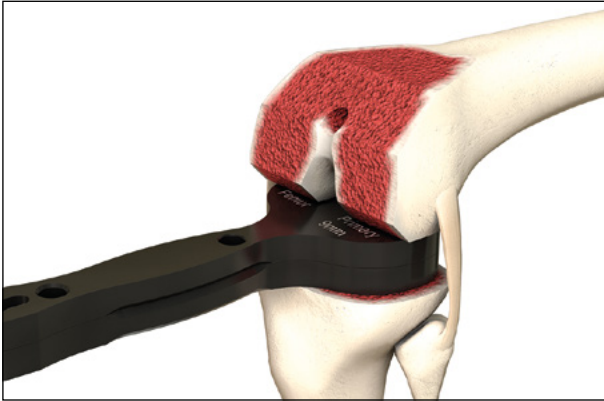


Fig. 32 Inspection of the flexion gap

Inspection of the flexion gap

- Insertion of the Spacer Block into the flexion gap with the Spacer Block Insert, previously defined in extension (Fig. 32)
- Evaluation of ligament stability, both medial and lateral

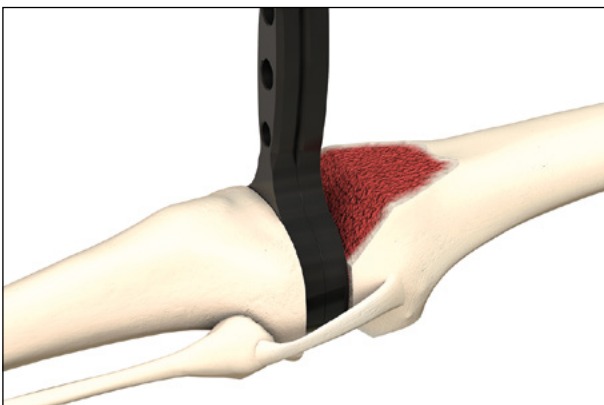


Fig. 33 Inspection of the extension gap

Inspection of the extension gap

Verification of the extension gap by inserting the Spacer Block with the appropriate Spacer Block Insert (Fig. 33).

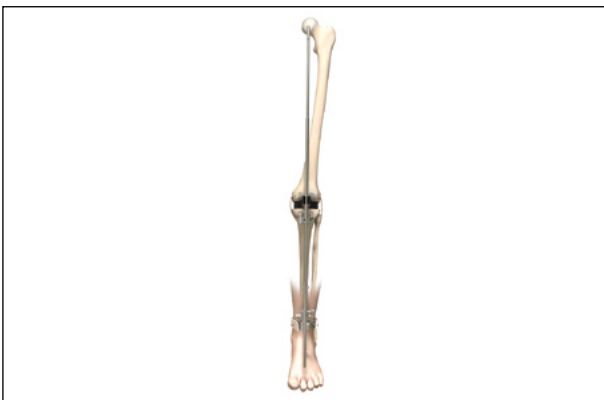


Fig. 34

Check the mechanical axis (Fig. 34), medial and lateral stability as well as extension ability. If conditions are too tight, a re-osteotomy can be performed at the distal femur or on the tibia.

Removal of the Spacer Block.

Surgical technique

3. Preparation and implantation

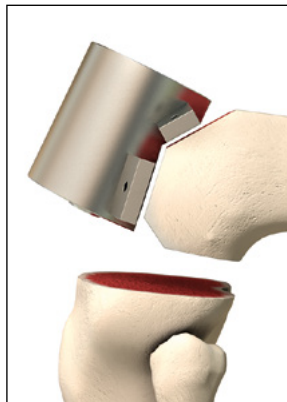


Fig. 35 Positioning of the Trochlea Bushing

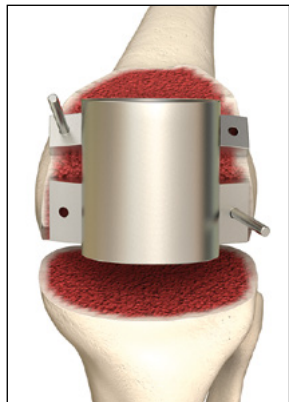


Fig. 36 Fixation of the Trochlea Bushing



Fig. 37

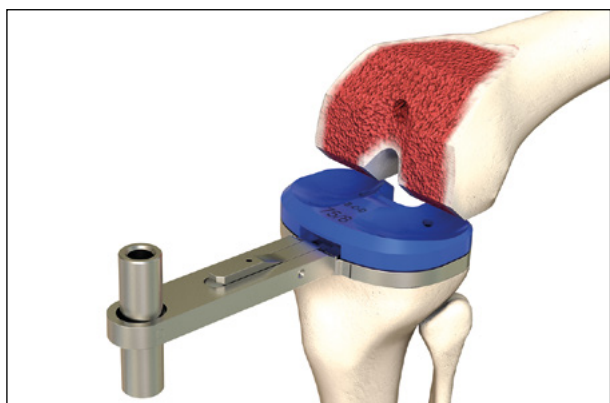


Fig. 38 Insertion of the Tibial Trial Prosthesis

3.1 balanSys CR, UC and RP

Femoral preparation

The Trochlea Bushing is placed onto the femur (Fig. 35) and affixed with at least two Pins diagonally (Fig. 36).

Remark

The authors recommend placing the Trochlea Bushing slightly lateral for optimized patella tracking.

The trochlea is chamfered by pushing the Trochlear Reamer to the stop (Fig. 37).

Remove all instruments.

Insertion of the determined Tibial Template and PE Trial Inlay (Fig. 38).

Remark

It must be ensured that the chosen template provides the desired tibial coverage and takes the chosen femur size into account.

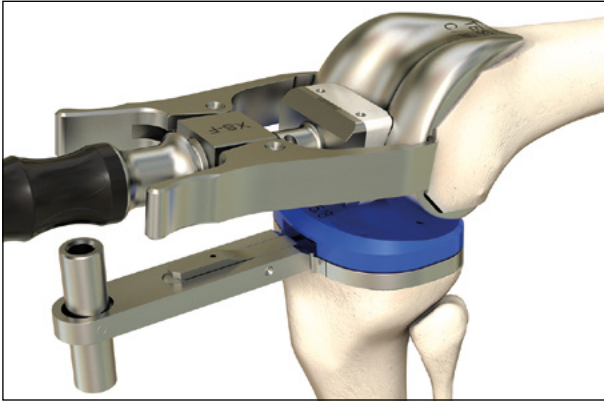


Fig. 39 Insertion of the Trial Femur

Insertion of the Trial Femur (Fig. 39) with the aid of the Femoral Holder and fixation with the Femoral Impactor (Fig. 40).

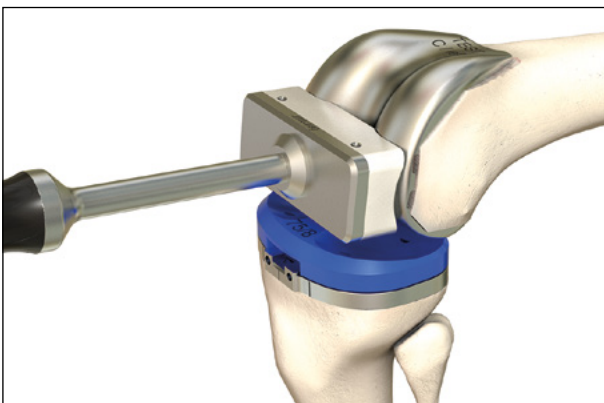


Fig. 40 Placement with the Femoral Impactor

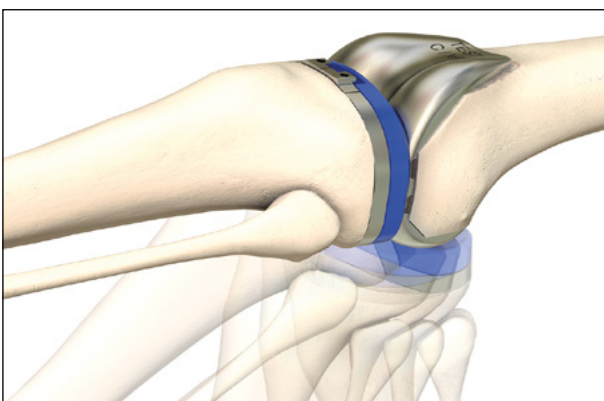


Fig. 41 Full movement test of the knee joint

Repositioning of the extension apparatus.

Remark

If the patella is replaced, it is recommended by the authors to perform the patellar osteotomy and to position the patellar trial component before the knee is tested.

With all the provisional components in place, the knee is tested with respect to maximal range of motion, stability, kinematics and mobility (Fig. 41).

Anterior marking of the position of the Tibial Template with the electrosurgical knife on the tibia.

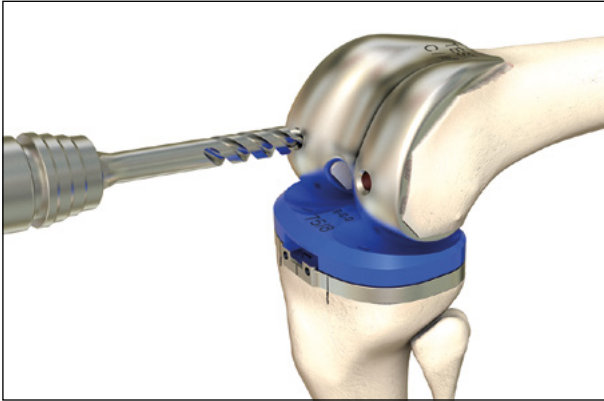


Fig. 42 Drilling for the femoral anchor pins

Drilling of the two holes for the femoral anchor pins (Fig. 42).

Removal of the Trial Femur and PE Trial Inlay.

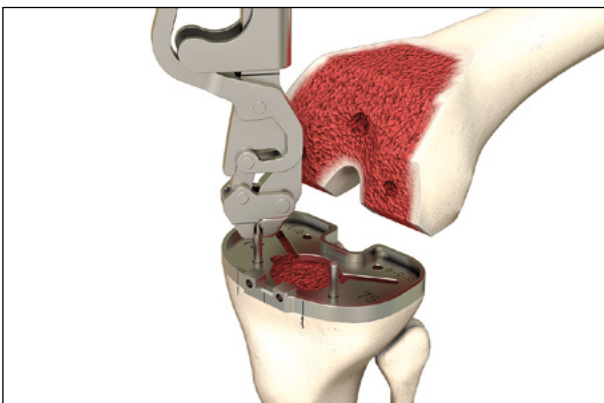


Fig. 43 Fixation of the Tibial Template

Tibial preparation

Fixation of the Tibial Template with two Pins (Fig. 43) matching the marking on the tibial head.



Fig. 44a

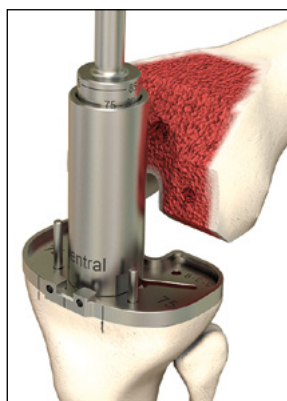


Fig. 44b

The Tibial Reamer must be connected to a power drill. Position the Reamer guide onto the Tibial Template. The three pegs need to be positioned in the openings of the Tibial Template. Note the marking **«ventral»**. Lock the Reamer Guide to the Tibial Template by turning the Reamer slightly clockwise until stop. In this position, the reamer Guide-Tibial Template assembly is locked (Fig. 44a). Drill out the tibial medullary space (Fig. 44b).

The depth must correspond to the appropriate length of the anchorage stem of the predetermined balanSys PS Tibial Plateau. Size markings are etched onto the Reamer.

Removal of the Reamer Guide.

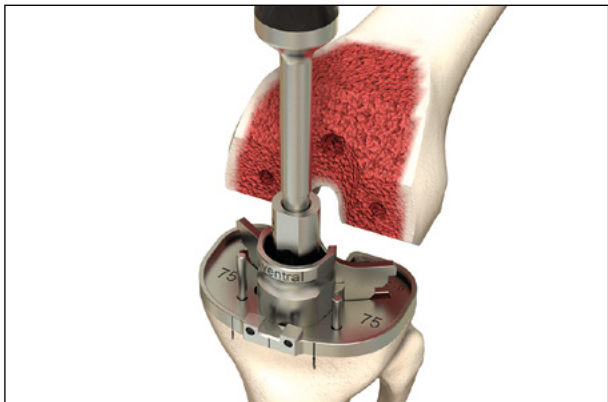


Fig. 45 Impaction of the Anchoring Bolt

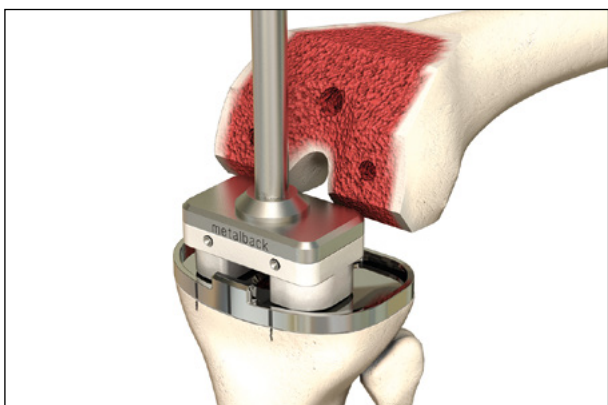


Fig. 46



Fig. 47

Assembling of the Fin Chisel with the Handle.

Position the Chisel Centring Guide onto the Tibial Template. The three pegs need to be positioned in the openings of the Tibial Template. Note the marking **«ventral»**. Lock the Chisel Centring Guide to the Tibial Template by turning it slightly clockwise until stop. In this position, the Chisel Centring Guide-Tibial Template assembly is locked. Impact the Fin Chisel assembly with care to prevent fracture on the tibia (Fig. 45).

During introduction of the Fin Chisel care has to be taken to protect the lateral collateral ligament and the popliteal tendon when using a surgical technique with a medial approach.

Impacting until the instrument bottoms out on the Tibial Template. The depths of the fins are defined by the size of the Tibial Template.

Removal of all remaining instruments.

Thorough cleaning of the osteotomy surfaces is required (e.g. with a pulse lavage).

Final implantation

The cement should be in the early dough phase when applied. Follow the instructions for the specific bone cement.

For secure fixation of the Tibial Plateau in the bone, it is necessary that the tibial backside be fully cemented in the dough phase of the cement. The stem and fins may be cemented or not.

Failure to fully cement and pressurize the Tibial Plateau may lead to early loosening of the prosthesis. Furthermore, cementing in advanced stages of polymerization can lead to early loosening of the prosthesis.

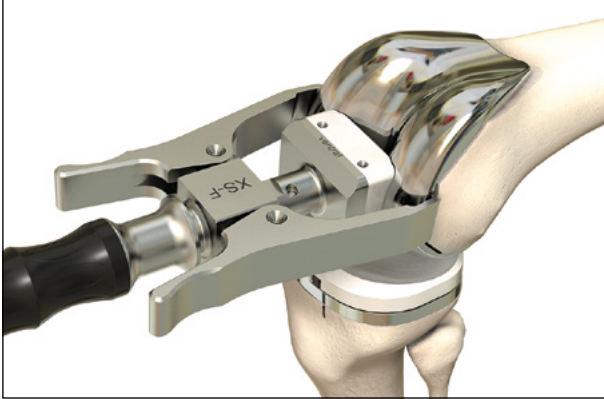


Fig. 48

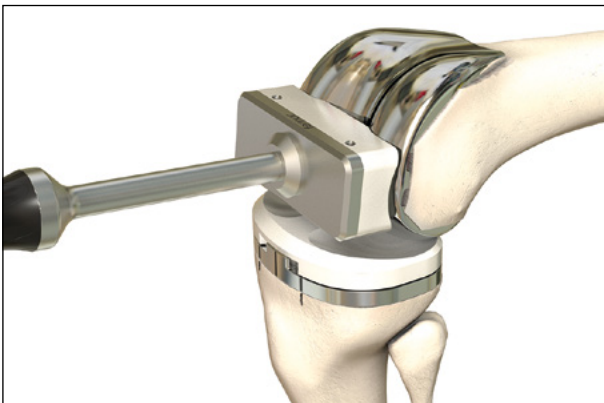


Fig. 49



Fig. 50

Remark

Careful application of cement should be considered to avoid excess cement on the posterior parts of the femur and femoral component. Cement that extrudes posteriorly is difficult to remove.

Insertion of the balanSys Tibial Plateau with the aid of the Tibial Impactor (Fig. 46).

Insertion of the balanSys Inlay (Fig. 47).

Insertion of the balanSys Femur (cemented or uncemented) with the aid of the Femoral Holder (Fig. 48) and final impactation with the Femoral Impactor (Fig. 49).

The leg should be in extension during hardening of the bone cement (Fig. 50).

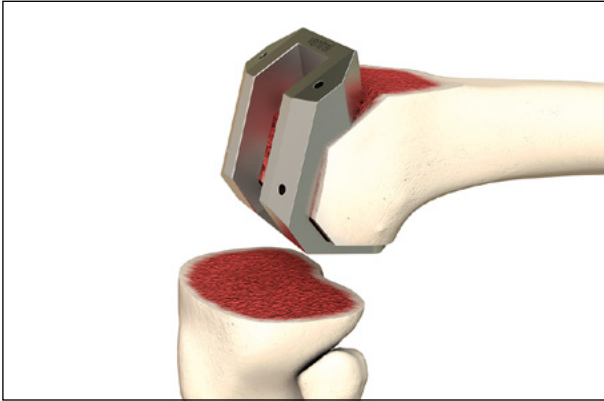


Fig. 51 Positioning of the Femur Box Cutting Guide

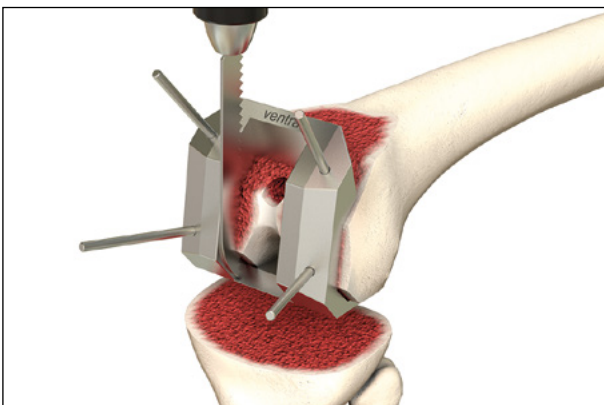


Fig. 52

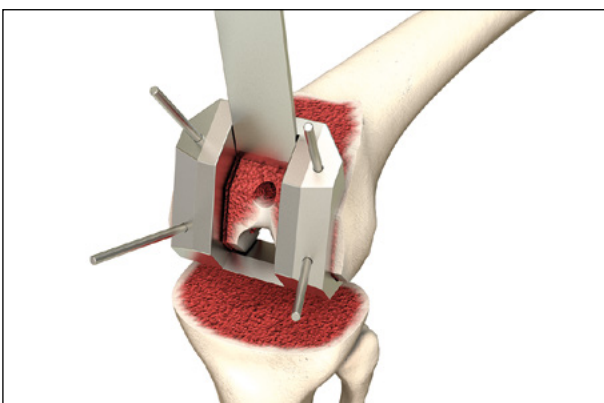


Fig. 53

3.2 balansys PS

Femoral preparation

The appropriate size of the Femur Box Cutting Guide is placed onto the femur. It must rest flush onto the resected surfaces of the posterior and distal femur (Fig. 51).

Remark

The authors recommend placing the Femur Box Cutting Guide slightly lateral for optimized patella tracking.

The guide will contact the resected bone surfaces and must be secured to the femur with four Pins diagonally. The most posterior Pins have to be introduced first.

A reciprocating saw should be used and guided along the walls of the open box to cut the medial and lateral sides and the base of the intercondylar notch (Fig. 52).

The base of the intercondylar notch can optionally be cut out with the balansys Chisel (Fig. 53).

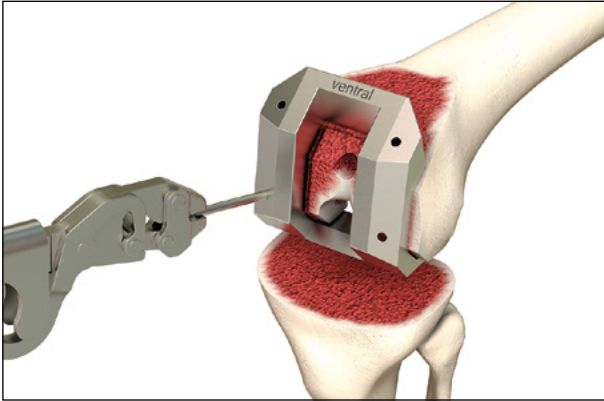


Fig. 54 Removal of Pins with Pliers

After the cuts, the Femur Box Cutting Guide and the Pins are removed (Fig. 54) in the following order:

1. Pins
2. Cutting guide
3. Resected bone block

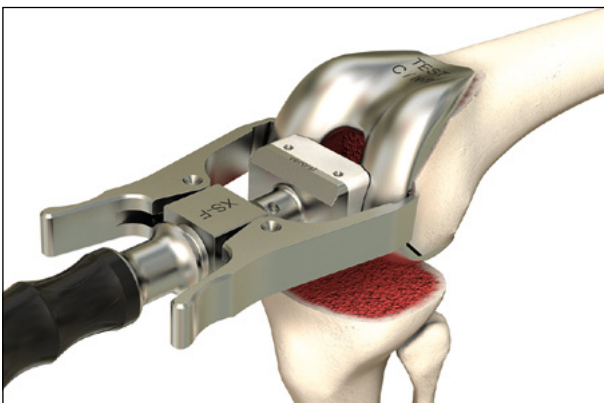


Fig. 55 Insertion of the Trial Femur

Determination of final implant position

Insertion of the Trial Femur with the Femoral Holder (Fig. 55).

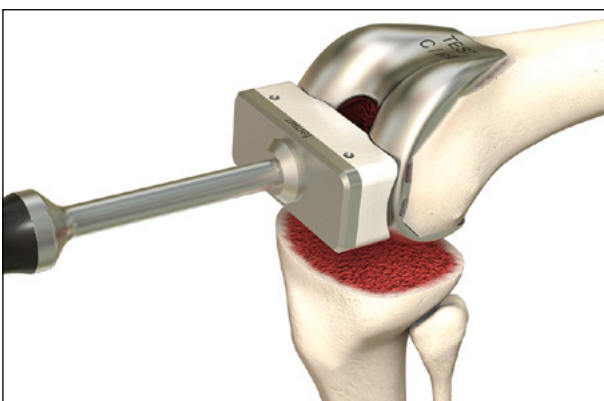


Fig. 56 Fixation with the Femoral Impactor

Final fixation must be performed with the Femoral Impactor (Fig. 56).

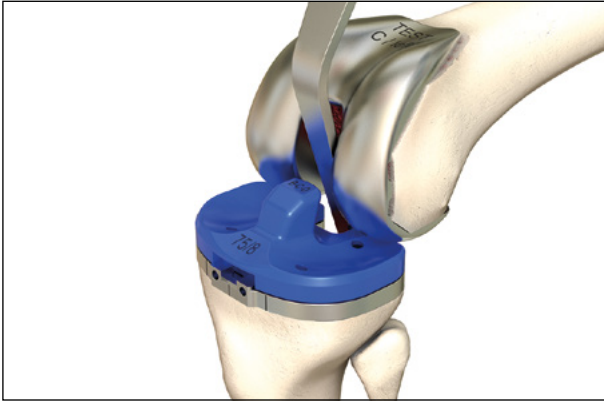


Fig. 57 Insertion of the Tibial Template and PE Trial Inlay

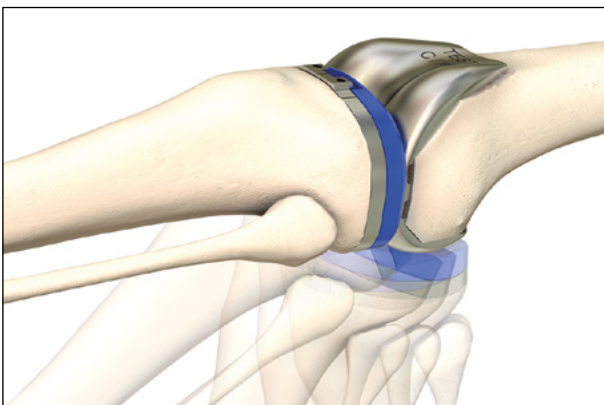


Fig. 58 Full movement test of the knee joint

Once the femur is impacted, the tibia has to be subluxed anteriorly with Hohmann retractors in order to position the Tibial Template and the PE Trial Inlay (Fig. 57).

Insertion of the Tibial Template and the PE Trial Inlay in the predetermined dimension and positioning onto the resected tibia.

Remark

It must be ensured that the chosen template provides the desired tibial coverage and takes the chosen femur size into account.

Repositioning of the extension apparatus.

Remark

If the patella is replaced, it is recommended by the authors to perform the patellar osteotomy and to position the patellar trial component before the knee is tested.

With all the provisional PS components in place, the knee is tested with respect to maximal range of motion, stability, kinematics and mobility (Fig. 58).

The tibial component tends to seat itself in the position where it best articulates with the femur. Any desired adjustments should be made before the final position is marked.

After final placement has been determined, a cautery is used to mark the location of the laser etch lines on the anterior periphery of the Tibial Template.

Removal of the PE Trial Inlay and the Trial Femur.



Fig. 59 Fixation with Pins

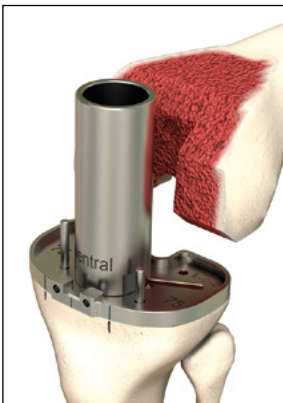


Fig. 60 Position of the Reamer Guide

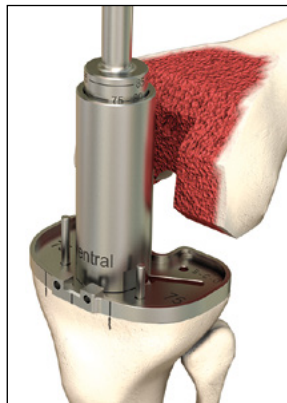


Fig. 61 Drill out of the tibial medullary space

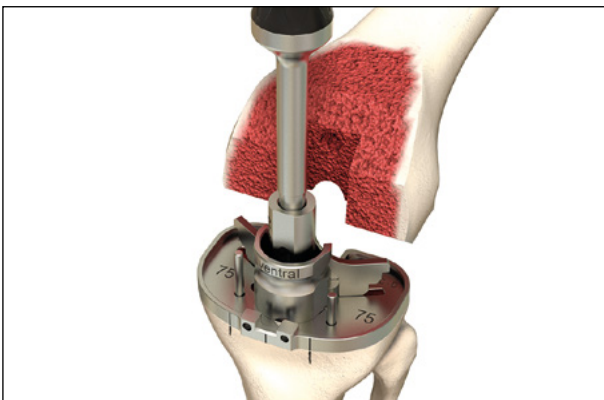


Fig. 62 Preparation of the fin

Tibial preparation

Re-positioning of the Tibial Template according to the marking of the cautery and fixation with two anterior Pins (Fig. 59).

The Tibial Reamer must be connected to a power drill.

Position the Reamer guide onto the Tibial Template. The three pegs need to be positioned in the openings of the Tibial Template. Note the marking **«ventral»**. Lock the Reamer Guide to the Tibial Template by turning the Reamer slightly clockwise until stop. In this position, the Reamer Guide-Tibial Template assembly is locked (Fig. 60). Drill out the tibial medullary space (Fig. 61).

The depth must correspond to the appropriate length of the anchorage stem of the predetermined balanSys PS Tibial Plateau. Size markings are etched onto the Reamer.

Removal of the Reamer Guide.

Assembling of the Fin Chisel with the Handle.

Position the Chisel Centring Guide onto the Tibial Template. The three pegs need to be positioned in the openings of the Tibial Template. Note the marking **«ventral»**. Lock the Chisel Centring Guide to the Tibial Template by turning it slightly clockwise until stop. In this position, the Chisel Centring Guide-Tibial Template assembly is locked. Impact the fin Chisel assembly with care to prevent fracture on the tibia (Fig. 62).

During introduction of the Fin Chisel care has to be taken to protect the lateral collateral ligament and the popliteal tendon when using a surgical technique with a medial approach.

Impacting until the instrument bottoms out on the Tibial Template. The depths of the fins are defined by the size of the Tibial Template.

Removal of all remaining instruments.

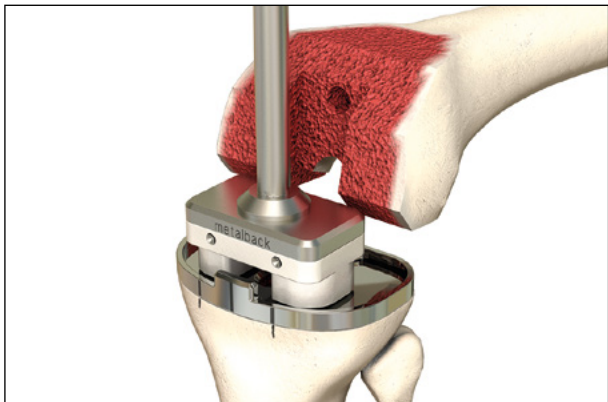


Fig. 63 Insert tibial component and impact it

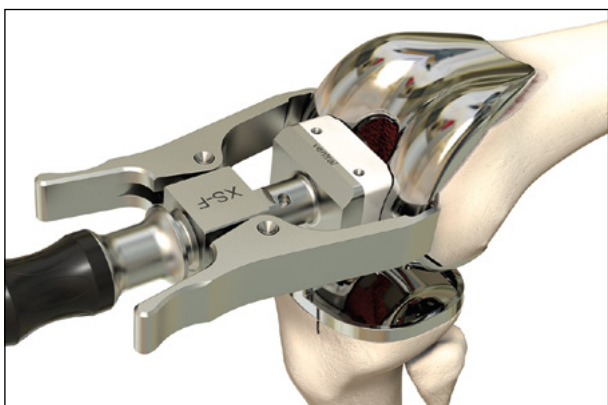


Fig. 64 Insertion of the femur

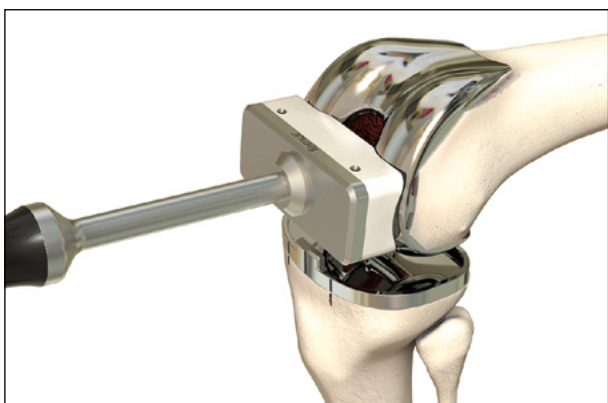


Fig. 65 Final impactation with the Femoral Impactor

Definitive implantation balanSys PS

After the implants have been chosen, one last check is recommended to ensure that the femoral, tibial, and inlay components match.

The cement should be in the early dough phase when applied. Follow the instructions for the specific bone cement.

For secure fixation of the Tibial Plateau in the bone, it is necessary that the tibial backside be fully cemented in the dough phase of the cement. The stem and fins may be cemented or not.

Failure to fully cement and pressurize the Tibial Plateau may lead to early loosening of the prosthesis. Furthermore, cementing in advanced stages of polymerization can lead to early loosening of the prosthesis.

Remark

Careful application of cement should be considered to avoid excess cement on the posterior parts of the femur and femoral component. Cement that extrudes posteriorly is difficult to remove.

For the final implantation, the knee should be in flexion.

Insert the balanSys PS Tibial Plateau and impact it with the aid of the Tibial Impactor (Fig. 63).

The balanSys PS Femur component is placed onto the femur with the aid of the Femoral Holder (Fig. 64). The femoral condyles of the balanSys PS Femur must be protected in order to prevent any scratching.

Final impactation with the Femoral Impactor (Fig. 65).

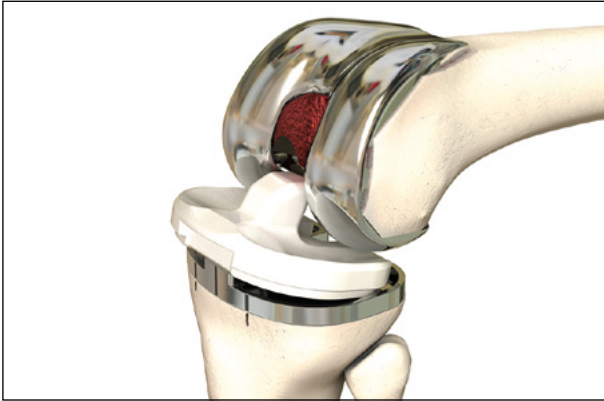


Fig. 66

Remark

Removal of excess cement. It is strongly recommended to give extra care to remove cement along the proximal portion of the femoral component and the femoral box. The surface of the tibial plateau must be cleaned.

Insertion of the balanSys PS Inlay (Fig. 66).



Fig. 67

The leg should be in extension during hardening of the bone cement (Fig. 67).

Remark

After insertion of the Inlay and hardening of the bone cement, extra care should be given to remove any remaining cement inside the femoral box.

Appendix

1 – Intramedullary tibial osteotomy

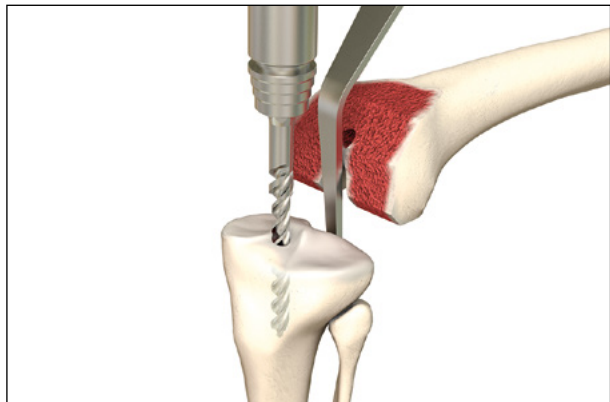


Fig. 68 Opening of the medullary canal

The tibial osteotomy can also be performed using intramedullary alignment.

The medullary canal is opened in the longitudinal axis at the eminentia with a 8.5 mm Drill (Fig. 68).

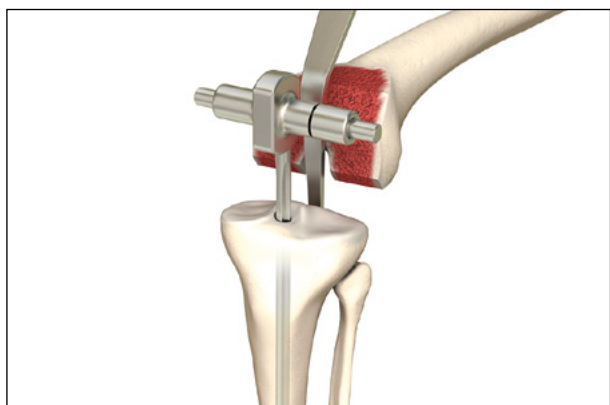


Fig. 69

Slow insertion of the Trs. Intramedullary Rod (Fig. 69).

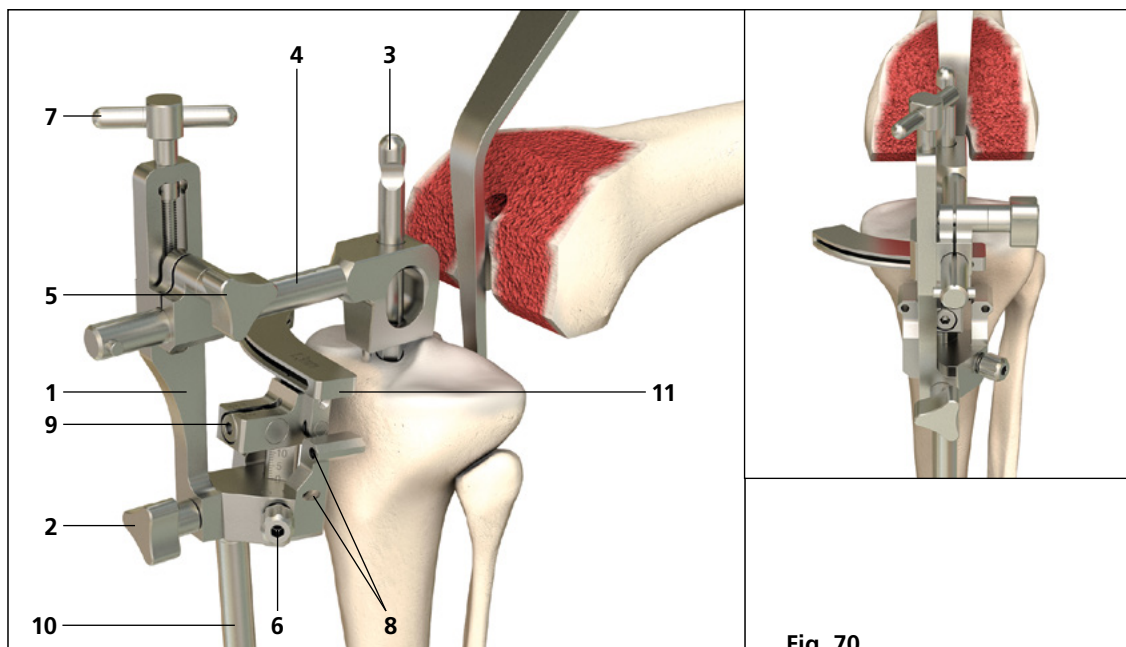


Fig. 70

Positioning of the intramedullary aligned Tibial Reference System (Fig. 70)

- Connect the intramedullary Tibial Reference System (1) with the Trs. Aiming Device Proximal (10)
- Fix it with the screw (2)
- Slide pre-mounted equipment onto the Trs. Intramedullary Rod (3)
- Distal alignment of the Tibial Reference System to the second metatarsophalangeal bone and proximally to the transition of the medial to the middle third of the tibial tuberosity
- Impaction on the bracket (4)
- Adjustment of the varus-valgus angle with the aid of the adjusting screw (5)
- First set the Tibial Cutting Guide (11) to zero and fix it with the screw (6). Determine the slope. Use the Hexagonal Screwdriver to set the posterior slope in accordance with the anatomy, as required from 0° – 11° (9)
- Determine the original joint line at the height of the best preserved tibial joint surface with the Reference Plate or Tibial Stylus
- The adjustment of the joint line is made with the aid of the T-screw (7)
- Fix the reference system proximally with at least two Pins in the specified holes (oblique and straight) (8)
- The holes are pre-drilled with the 3.2 mm Drill.
- Adjust the Tibial Cutting Guide 6–8mm distally with the aid of the Hexagonal Screwdriver (6)
- Check the planned tibial osteotomy with the aid of the Reference Plate
- Removal screw of the intramedullary retainer (2) and the intramedullary aligned tibial resection instrument (1, 3, 4)

Remark

Our authors recommend a posterior slope of 7° for a cruciate-retaining implant and up to 5° for a PCL-substituting implant.

For continuation see section «Tibial osteotomy», (Fig. 29).

Appendix

2 – balanSys 3-Peg Patella

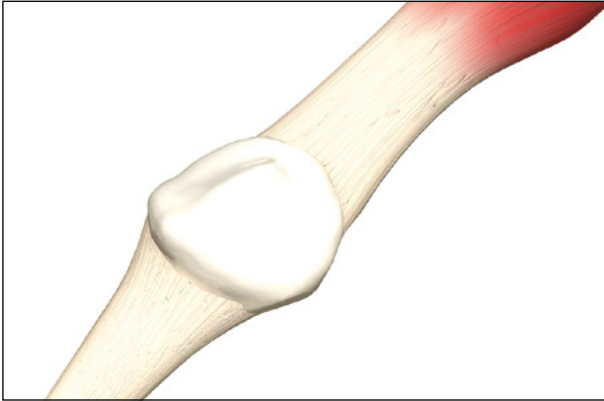


Fig. 71

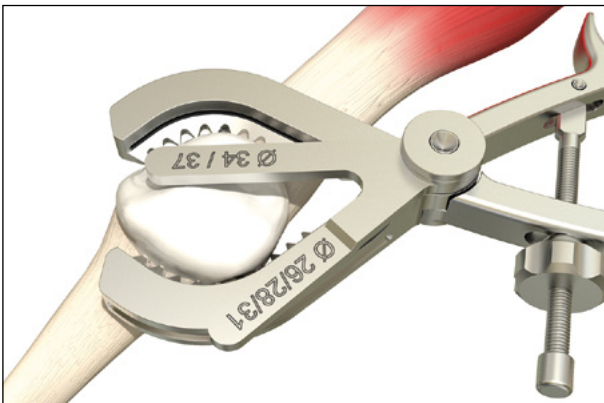


Fig. 72

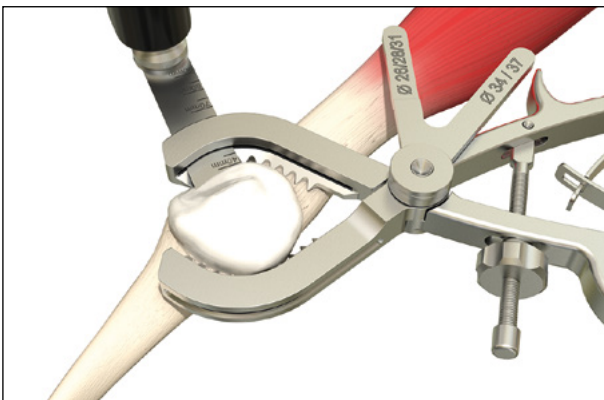


Fig. 73

Remove any marginal osteophytes.

Determine the size of the patella on the basis of the anatomical situation.

Dia	3-Peg FLAT	3-Peg
26	8 mm	–
28	8 mm	10.2 mm
31	8 mm	11.4 mm
34	9 mm	12.3 mm
37	9 mm	13.0 mm

Evert the patella and grip centrally with the Patella Forceps.

Adjust the resection height with the height topper to the selected patellar size setting.

Remark

Ensure that you are using the Patella Resection Pliers Flat for the 3-Peg Patella FLAT (marked with sizes **26**–**37**) and the Patella Resection Pliers Standard for the 3-Peg Patella (marked with sizes **28**–**37**), respectively.

Perform the patellar osteotomy through the saw guide on lateral side of the Patella Forceps.

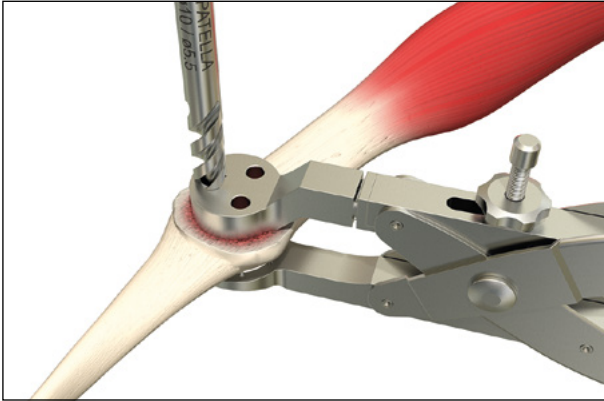


Fig. 74

Place the Drill Guide in order to determine the final positioning of the patella implant with regard to the predetermined gliding path of the femoral shield.

Drill the three patella peg holes with the 5.5 mm Drill Bit.



Fig. 75

Retro patellar surface prepared for implantation.



Fig. 76

Insert the Trial Patella of predetermined size (FLAT or standard).

Chamfer the medial and lateral edges of the rear surface of the patella.

Check the gliding movement in the femoro-patellar joint with respect to centering.

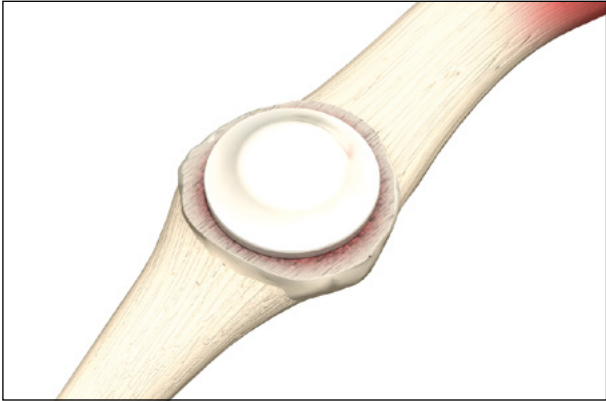


Fig. 77

Thoroughly clean the osteotomy surface.

Insertion of balanSys 3-Peg Patella cemented using the cement clamp.

Remove any marginal cement rests.

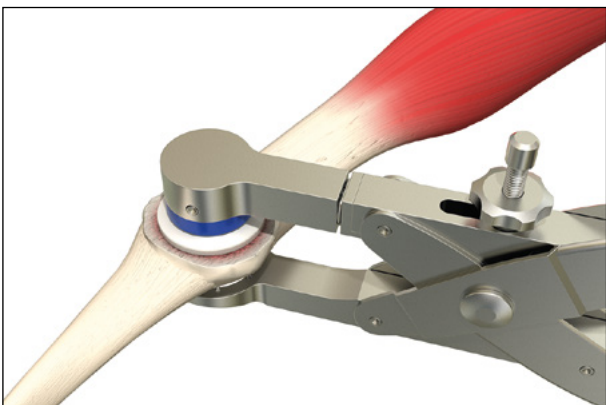


Fig. 78

After hardening of the cement remove the Cement Clamp.

After closing the joint capsule, make a final functional test and check the centred running of the patella.

Appendix

3 – Item numbers of the balanSys implants

balanSys Femur Components for CR/UC/RP

balanSys Femur, cemented



Item no.	Mediolat.	Size
72.15.3401	56mm	XS left
72.15.3701	58mm	S left
72.15.4001	60mm	A left
72.15.4301	64mm	B left
72.15.4601	68mm	C left
72.15.4901	72mm	D left
72.15.5201	76mm	E left
72.15.5501	80mm	F left
72.15.3402	56mm	XS right
72.15.3702	58mm	S right
72.15.4002	60mm	A right
72.15.4302	64mm	B right
72.15.4602	68mm	C right
72.15.4902	72mm	D right
72.15.5202	76mm	E right
72.15.5502	80mm	F right

Material: CoCrMo

balanSys Femur, uncemented



Item no.	Mediolat.	Size
73.15.3401TPS	56mm	XS left
73.15.3701TPS	58mm	S left
73.15.4001TPS	60mm	A left
73.15.4301TPS	64mm	B left
73.15.4601TPS	68mm	C left
73.15.4901TPS	72mm	D left
73.15.5201TPS	76mm	E left
73.15.5501TPS	80mm	F left
73.15.3402TPS	56mm	XS right
73.15.3702TPS	58mm	S right
73.15.4002TPS	60mm	A right
73.15.4302TPS	64mm	B right
73.15.4602TPS	68mm	C right
73.15.4902TPS	72mm	D right
73.15.5202TPS	76mm	E right
73.15.5502TPS	80mm	F right

Material: CoCrMo, TiCP coated

balanSys Fixed Bearing Components

balanSys CR PE Inlay



Item no.	Mediolat.	Size
74.30.5908	59mm	8.0mm
74.30.5910	59mm	10.5mm
74.30.5913	59mm	13.0mm
74.30.5915	59mm	15.5mm
74.30.6208	62mm	8.0mm
74.30.6210	62mm	10.5mm
74.30.6213	62mm	13.0mm
74.30.6215	62mm	15.5mm
74.30.6408	64mm	8.0mm
74.30.6410	64mm	10.5mm
74.30.6413	64mm	13.0mm
74.30.6415	64mm	15.5mm
74.30.6708	67mm	8.0mm
74.30.6710	67mm	10.5mm
74.30.6713	67mm	13.0mm
74.30.6715	67mm	15.5mm

Item no.	Mediolat.	Size
74.30.7008	70mm	8.0mm
74.30.7010	70mm	10.5mm
74.30.7013	70mm	13.0mm
74.30.7015	70mm	15.5mm
74.30.7508	75mm	8.0mm
74.30.7510	75mm	10.5mm
74.30.7513	75mm	13.0mm
74.30.7515	75mm	15.5mm
72.34.0170	80mm	8.0mm
72.34.0171	80mm	10.5mm
72.34.0172	80mm	13.0mm
72.34.0173	80mm	15.5mm
72.34.0174	85mm	8.0mm
72.34.0175	85mm	10.5mm
72.34.0176	85mm	13.0mm
72.34.0177	85mm	15.5mm

Material: UHMWPE



balanSys CR vitamys Inlay

Item no.	Mediolat.	Size
72.34.1000	59mm	8.0mm
72.34.1001	59mm	9.0mm
72.34.1002	59mm	10.5mm
72.34.1003	59mm	11.5mm
72.34.1004	59mm	13.0mm
72.34.1005	59mm	15.5mm
72.34.1010	62mm	8.0mm
72.34.1011	62mm	9.0mm
72.34.1012	62mm	10.5mm
72.34.1013	62mm	11.5mm
72.34.1014	62mm	13.0mm
72.34.1015	62mm	15.5mm
72.34.1020	64mm	8.0mm
72.34.1021	64mm	9.0mm
72.34.1022	64mm	10.5mm
72.34.1023	64mm	11.5mm
72.34.1024	64mm	13.0mm
72.34.1025	64mm	15.5mm
72.34.1030	67mm	8.0mm
72.34.1031	67mm	9.0mm
72.34.1032	67mm	10.5mm
72.34.1033	67mm	11.5mm
72.34.1034	67mm	13.0mm
72.34.1035	67mm	15.5mm

Material: VEPE

Item no.	Mediolat.	Size
72.34.1040	70mm	8.0mm
72.34.1041	70mm	9.0mm
72.34.1042	70mm	10.5mm
72.34.1043	70mm	11.5mm
72.34.1044	70mm	13.0mm
72.34.1045	70mm	15.5mm
72.34.1050	75mm	8.0mm
72.34.1051	75mm	9.0mm
72.34.1052	75mm	10.5mm
72.34.1053	75mm	11.5mm
72.34.1054	75mm	13.0mm
72.34.1055	75mm	15.5mm
72.34.1060	80mm	8.0mm
72.34.1061	80mm	9.0mm
72.34.1062	80mm	10.5mm
72.34.1063	80mm	11.5mm
72.34.1064	80mm	13.0mm
72.34.1065	80mm	15.5mm
72.34.1070	85mm	8.0mm
72.34.1071	85mm	9.0mm
72.34.1072	85mm	10.5mm
72.34.1073	85mm	11.5mm
72.34.1074	85mm	13.0mm
72.34.1075	85mm	15.5mm



balanSys UC PE Inlay

Item no.	Mediolat.	Size
77.30.5908	59mm	8.0mm
77.30.5910	59mm	10.5mm
77.30.5913	59mm	13.0mm
77.30.5915	59mm	15.5mm
77.30.5918	59mm	18.0mm
77.30.5920	59mm	20.5mm
77.30.6208	62mm	8.0mm
77.30.6210	62mm	10.5mm
77.30.6213	62mm	13.0mm
77.30.6215	62mm	15.5mm
77.30.6218	62mm	18.0mm
77.30.6220	62mm	20.5mm
77.30.6408	64mm	8.0mm
77.30.6410	64mm	10.5mm
77.30.6413	64mm	13.0mm
77.30.6415	64mm	15.5mm
77.30.6418	64mm	18.0mm
77.30.6420	64mm	20.5mm
77.30.6708	67mm	8.0mm
77.30.6710	67mm	10.5mm
77.30.6713	67mm	13.0mm
77.30.6715	67mm	15.5mm
77.30.6718	67mm	18.0mm
77.30.6720	67mm	20.5mm

Material: UHMWPE

Item no.	Mediolat.	Size
77.30.7008	70mm	8.0mm
77.30.7010	70mm	10.5mm
77.30.7013	70mm	13.0mm
77.30.7015	70mm	15.5mm
77.30.7018	70mm	18.0mm
77.30.7020	70mm	20.5mm
77.30.7508	75mm	8.0mm
77.30.7510	75mm	10.5mm
77.30.7513	75mm	13.0mm
77.30.7515	75mm	15.5mm
77.30.7518	75mm	18.0mm
77.30.7520	75mm	20.5mm
72.34.0182	80mm	8.0mm
72.34.0183	80mm	10.5mm
72.34.0184	80mm	13.0mm
72.34.0185	80mm	15.5mm
72.34.0186	80mm	18.0mm
72.34.0187	80mm	20.5mm
72.34.0188	85mm	8.0mm
72.34.0189	85mm	10.5mm
72.34.0190	85mm	13.0mm
72.34.0191	85mm	15.5mm
72.34.0192	85mm	18.0mm
72.34.0193	85mm	20.5mm



vitamys®

balanSys UC vitamys Inlay

Item no.	Mediolat.	Size
72.34.1100	59mm	8.0mm
72.34.1101	59mm	9.0mm
72.34.1102	59mm	10.5mm
72.34.1103	59mm	11.5mm
72.34.1104	59mm	13.0mm
72.34.1105	59mm	15.5mm
72.34.1106	59mm	18.0mm
72.34.1107*	59mm	20.5mm
72.34.1110	62mm	8.0mm
72.34.1111	62mm	9.0mm
72.34.1112	62mm	10.5mm
72.34.1113	62mm	11.5mm
72.34.1114	62mm	13.0mm
72.34.1115	62mm	15.5mm
72.34.1116	62mm	18.0mm
72.34.1117*	62mm	20.5mm
72.34.1120	64mm	8.0mm
72.34.1121	64mm	9.0mm
72.34.1122	64mm	10.5mm
72.34.1123	64mm	11.5mm
72.34.1124	64mm	13.0mm
72.34.1125	64mm	15.5mm
72.34.1126	64mm	18.0mm
72.34.1127*	64mm	20.5mm
72.34.1130	67mm	8.0mm
72.34.1131	67mm	9.0mm
72.34.1132	67mm	10.5mm
72.34.1133	67mm	11.5mm
72.34.1134	67mm	13.0mm
72.34.1135	67mm	15.5mm
72.34.1136	67mm	18.0mm
72.34.1137*	67mm	20.5mm

Material: VEPE

*on request

Item no.	Mediolat.	Size
72.34.1140	70mm	8.0mm
72.34.1141	70mm	9.0mm
72.34.1142	70mm	10.5mm
72.34.1143	70mm	11.5mm
72.34.1144	70mm	13.0mm
72.34.1145	70mm	15.5mm
72.34.1146	70mm	18.0mm
72.34.1147*	70mm	20.5mm
72.34.1150	75mm	8.0mm
72.34.1151	75mm	9.0mm
72.34.1152	75mm	10.5mm
72.34.1153	75mm	11.5mm
72.34.1154	75mm	13.0mm
72.34.1155	75mm	15.5mm
72.34.1156	75mm	18.0mm
72.34.1157*	75mm	20.5mm
72.34.1160	80mm	8.0mm
72.34.1161	80mm	9.0mm
72.34.1162	80mm	10.5mm
72.34.1163	80mm	11.5mm
72.34.1164	80mm	13.0mm
72.34.1165	80mm	15.5mm
72.34.1166	80mm	18.0mm
72.34.1167*	80mm	20.5mm
72.34.1170	85mm	8.0mm
72.34.1171	85mm	9.0mm
72.34.1172	85mm	10.5mm
72.34.1173	85mm	11.5mm
72.34.1174	85mm	13.0mm
72.34.1175	85mm	15.5mm
72.34.1176	85mm	18.0mm
72.34.1177*	85mm	20.5mm



balanSys PS Tibial Plateau, cemented

Item no.	Mediolateral
79.15.0400	59 mm
79.15.0401	62 mm
79.15.0056	64 mm
79.15.0402	67 mm
79.15.0057	70 mm
79.15.0058	75 mm
79.15.0059	80 mm
79.15.0060	85 mm

Material: CoCrMo

balanSys Mobile Bearing RP Components

balanSys RP PE Inlay



Item no.	Femur	Size	Item no.	Femur	Size
72.34.0200	XS	8.0 mm	78.30.7008	C	8.0 mm
72.34.0201	XS	10.5 mm	78.30.7010	C	10.5 mm
72.34.0202	XS	13.0 mm	78.30.7013	C	13.0 mm
72.34.0203	XS	15.5 mm	78.30.7015	C	15.5 mm
72.34.0206	S	8.0 mm	78.30.7408	D	8.0 mm
72.34.0207	S	10.5 mm	78.30.7410	D	10.5 mm
72.34.0208	S	13.0 mm	78.30.7413	D	13.0 mm
72.34.0209	S	15.5 mm	78.30.7415	D	15.5 mm
78.30.6208	A	8.0 mm	78.30.7808	E	8.0 mm
78.30.6210	A	10.5 mm	78.30.7810	E	10.5 mm
78.30.6213	A	13.0 mm	78.30.7813	E	13.0 mm
78.30.6215	A	15.5 mm	78.30.7815	E	15.5 mm
78.30.6608	B	8.0 mm	72.34.0242	F	8.0 mm
78.30.6610	B	10.5 mm	72.34.0243	F	10.5 mm
78.30.6613	B	13.0 mm	72.34.0244	F	13.0 mm
78.30.6615	B	15.5 mm	72.34.0245	F	15.5 mm

Material: UHMWPE, FeCrNiMoMn (Contrast balls, optional)



vitamys®

balanSys RP vitamys Inlay

Item no.	Femur	Size
72.34.1200	XS	8.0 mm
72.34.1201	XS	9.0 mm
72.34.1202	XS	10.5 mm
72.34.1203	XS	11.5 mm
72.34.1204	XS	13.0 mm
72.34.1205	XS	15.5 mm
72.34.1210	S	8.0 mm
72.34.1211	S	9.0 mm
72.34.1212	S	10.5 mm
72.34.1213	S	11.5 mm
72.34.1214	S	13.0 mm
72.34.1215	S	15.5 mm
72.34.1220	A	8.0 mm
72.34.1221	A	9.0 mm
72.34.1222	A	10.5 mm
72.34.1223	A	11.5 mm
72.34.1224	A	13.0 mm
72.34.1225	A	15.5 mm
72.34.1230	B	8.0 mm
72.34.1231	B	9.0 mm
72.34.1232	B	10.5 mm
72.34.1233	B	11.5 mm
72.34.1234	B	13.0 mm
72.34.1235	B	15.5 mm

Material: VEPE

Item no.	Femur	Size
72.34.1240	C	8.0 mm
72.34.1241	C	9.0 mm
72.34.1242	C	10.5 mm
72.34.1243	C	11.5 mm
72.34.1244	C	13.0 mm
72.34.1245	C	15.5 mm
72.34.1250	D	8.0 mm
72.34.1251	D	9.0 mm
72.34.1252	D	10.5 mm
72.34.1253	D	11.5 mm
72.34.1254	D	13.0 mm
72.34.1255	D	15.5 mm
72.34.1260	E	8.0 mm
72.34.1261	E	9.0 mm
72.34.1262	E	10.5 mm
72.34.1263	E	11.5 mm
72.34.1264	E	13.0 mm
72.34.1265	E	15.5 mm
72.34.1270	F	8.0 mm
72.34.1271	F	9.0 mm
72.34.1272	F	10.5 mm
72.34.1273	F	11.5 mm
72.34.1274	F	13.0 mm
72.34.1275	F	15.5 mm



balanSys RP Tibial Plateau, cemented

Item no.	Mediolateral
72.34.0059	59 mm
72.34.0060	62 mm
72.34.0061	64 mm
72.34.0062	67 mm
72.34.0063	70 mm
72.34.0064	75 mm
72.34.0065	80 mm
72.34.0066	85 mm

Material: CoCrMo



balanSys PS Components

balanSys PS Femur, cemented

Item no.	Mediolat.	Size
79.15.0999	56mm	XS right
79.15.1000	58mm	S right
79.15.0001	60mm	A right
79.15.0002	64mm	B right
79.15.0003	68mm	C right
79.15.0004	72mm	D right
79.15.0005	76mm	E right
79.15.1006	80mm	F right
79.15.1009	56mm	XS left
79.15.1010	58mm	S left
79.15.0011	60mm	A left
79.15.0012	64mm	B left
79.15.0013	68mm	C left
79.15.0014	72mm	D left
79.15.0015	76mm	E left
79.15.1016	80mm	F left

Material: CoCrMo



balanSys PS PE Inlay

Item no.	Mediolat.	Size
79.30.9986	59mm	8.0mm
79.30.9987	59mm	10.5mm
79.30.9988	59mm	13.0mm
79.30.9989	59mm	15.5mm
79.30.9990	59mm	18.0mm
79.30.9991	59mm	20.5mm
79.30.9992	59mm	23.0mm
79.30.9993	62mm	8.0mm
79.30.9994	62mm	10.5mm
79.30.9995	62mm	13.0mm
79.30.9996	62mm	15.5mm
79.30.9997	62mm	18.0mm
79.30.9998	62mm	20.5mm
79.30.9999	62mm	23.0mm
79.30.0200	64mm	8.0mm
79.30.0201	64mm	10.5mm
79.30.0202	64mm	13.0mm
79.30.0203	64mm	15.5mm
79.30.0204	64mm	18.0mm
79.30.0205	64mm	20.5mm
79.30.0206	64mm	23.0mm
79.30.0210	67mm	8.0mm
79.30.0211	67mm	10.5mm
79.30.0212	67mm	13.0mm
79.30.0213	67mm	15.5mm
79.30.0214	67mm	18.0mm
79.30.0215	67mm	20.5mm
79.30.0216	67mm	23.0mm

Material: UHMWPE

Item no.	Mediolat.	Size
79.30.0010	70mm	8.0mm
79.30.0011	70mm	10.5mm
79.30.0012	70mm	13.0mm
79.30.0013	70mm	15.5mm
79.30.0014	70mm	18.0mm
79.30.0015	70mm	20.5mm
79.30.0016	70mm	23.0mm
79.30.0020	75mm	8.0mm
79.30.0021	75mm	10.5mm
79.30.0022	75mm	13.0mm
79.30.0023	75mm	15.5mm
79.30.0024	75mm	18.0mm
79.30.0025	75mm	20.5mm
79.30.0026	75mm	23.0mm
72.34.0255	80mm	8.0mm
72.34.0256	80mm	10.5mm
72.34.0257	80mm	13.0mm
72.34.0258	80mm	15.5mm
72.34.0259	80mm	18.0mm
72.34.0260	80mm	20.5mm
72.34.0261	80mm	23.0mm
72.34.0262	85mm	8.0mm
72.34.0263	85mm	10.5mm
72.34.0264	85mm	13.0mm
72.34.0265	85mm	15.5mm
72.34.0266	85mm	18.0mm
72.34.0267	85mm	20.5mm
72.34.0268	85mm	23.0mm



balanSys PS vitamys Inlay

Item no.	Mediolat.	Size
72.34.1300	59mm	8.0mm
72.34.1301	59mm	9.0mm
72.34.1302	59mm	10.5mm
72.34.1303	59mm	11.5mm
72.34.1304	59mm	13.0mm
72.34.1305	59mm	15.5mm
72.34.1306	59mm	18.0mm
72.34.1307*	59mm	20.5mm
72.34.1310	62mm	8.0mm
72.34.1311	62mm	9.0mm
72.34.1312	62mm	10.5mm
72.34.1313	62mm	11.5mm
72.34.1314	62mm	13.0mm
72.34.1315	62mm	15.5mm
72.34.1316	62mm	18.0mm
72.34.1317*	62mm	20.5mm
72.34.1320	64mm	8.0mm
72.34.1321	64mm	9.0mm
72.34.1322	64mm	10.5mm
72.34.1323	64mm	11.5mm
72.34.1324	64mm	13.0mm
72.34.1325	64mm	15.5mm
72.34.1326	64mm	18.0mm
72.34.1327*	64mm	20.5mm
72.34.1330	67mm	8.0mm
72.34.1331	67mm	9.0mm
72.34.1332	67mm	10.5mm
72.34.1333	67mm	11.5mm
72.34.1334	67mm	13.0mm
72.34.1335	67mm	15.5mm
72.34.1336	67mm	18.0mm
72.34.1337*	67mm	20.5mm

Material: VEPE

*on request

Item no.	Mediolat.	Size
72.34.1340	70mm	8.0mm
72.34.1341	70mm	9.0mm
72.34.1342	70mm	10.5mm
72.34.1343	70mm	11.5mm
72.34.1344	70mm	13.0mm
72.34.1345	70mm	15.5mm
72.34.1346	70mm	18.0mm
72.34.1347*	70mm	20.5mm
72.34.1350	75mm	8.0mm
72.34.1351	75mm	9.0mm
72.34.1352	75mm	10.5mm
72.34.1353	75mm	11.5mm
72.34.1354	75mm	13.0mm
72.34.1355	75mm	15.5mm
72.34.1356	75mm	18.0mm
72.34.1357*	75mm	20.5mm
72.34.1360	80mm	8.0mm
72.34.1361	80mm	9.0mm
72.34.1362	80mm	10.5mm
72.34.1363	80mm	11.5mm
72.34.1364	80mm	13.0mm
72.34.1365	80mm	15.5mm
72.34.1366	80mm	18.0mm
72.34.1367*	80mm	20.5mm
72.34.1370	85mm	8.0mm
72.34.1371	85mm	9.0mm
72.34.1372	85mm	10.5mm
72.34.1373	85mm	11.5mm
72.34.1374	85mm	13.0mm
72.34.1375	85mm	15.5mm
72.34.1376	85mm	18.0mm
72.34.1377*	85mm	20.5mm

balanSys TiNbN Components

balanSys Femur TiNbN, cemented



Item no.	Mediolat.	Size
72.23.3401	56 mm	XS left
72.23.3701	58 mm	S left
72.23.4001	60 mm	A left
72.23.4301	64 mm	B left
72.23.4601	68 mm	C left
72.23.4901	72 mm	D left
72.23.5201	76 mm	E left
72.23.5501	80 mm	F left

Material: CoCrMo, TiNbN coating

Item no.	Mediolat.	Size
72.23.3402	56 mm	XS right
72.23.3702	58 mm	S right
72.23.4002	60 mm	A right
72.23.4302	64 mm	B right
72.23.4602	68 mm	C right
72.23.4902	72 mm	D right
72.23.5202	76 mm	E right
72.23.5502	80 mm	F right

balanSys PS Tibial Plateau TiNbN Fix, cemented



Item no.	Mediolateral
79.23.0400	59 mm
79.23.0401	62 mm
79.23.0056	64 mm
79.23.0402	67 mm

Material: CoCrMo, TiNbN coating

Item no.	Mediolateral
79.23.0057	70 mm
79.23.0058	75 mm
79.23.0059	80 mm
79.23.0060	85 mm

balanSys PS Femur TiNbN, cemented



Item no.	Mediolat.	Size
79.23.1009	56 mm	XS left
79.23.1010	58 mm	S left
79.23.0011	60 mm	A left
79.23.0012	64 mm	B left
79.23.0013	68 mm	C left
79.23.0014	72 mm	D left
79.23.0015	76 mm	E left
79.23.1016	80 mm	F left

Material: CoCrMo, TiNbN coating

Item no.	Mediolat.	Size
79.23.0999	56 mm	XS right
79.23.1000	58 mm	S right
79.23.0001	60 mm	A right
79.23.0002	64 mm	B right
79.23.0003	68 mm	C right
79.23.0004	72 mm	D right
79.23.0005	76 mm	E right
79.23.1006	80 mm	F right



balanSys 3-Peg Patella FLAT Components

Item no.	Diameter Ø
72.34.0049	26mm
72.34.0050	28mm
72.34.0051	31 mm
72.34.0052	34mm
72.34.0053	37mm

Material: UHMWPE, FeCrNiMoMn (Contrast balls)



balanSys 3-Peg Patella Components

Item no.	Diameter Ø
72.30.0128	28mm
72.30.0131	31 mm
72.30.0134	34mm
72.30.0137	37mm

Material: UHMWPE, FeCrNiMoMn (Contrast balls)

Appendix

5 – Item numbers of the balanSys instruments

Basic instruments

balanSys Basic Set 71.34.9000A	53
balanSys Basic Trial Set Essential 71.34.9005A	58

Surgical technique

balanSys LIS 4in1 71.34.9030A	61
-------------------------------	----

Trial instruments

balanSys Trial Set CR/UC	63
balanSys Trial Set PS	68
balanSys Trial Set RP (optional)	74


Patella instruments

balanSys 3-Peg Patella FLAT 71.34.0080A	78
balanSys 3-Peg Patella STANDARD 71.34.0081A	78

**Before each surgery, instruments should be checked for damage or deformation.
Use only undamaged instruments. Do not use trial components with marks or scratches.**

balanSys Basic Set 71.34.9000A

No Picture / 71.34.0545 **balanSys Basic Set Lid**

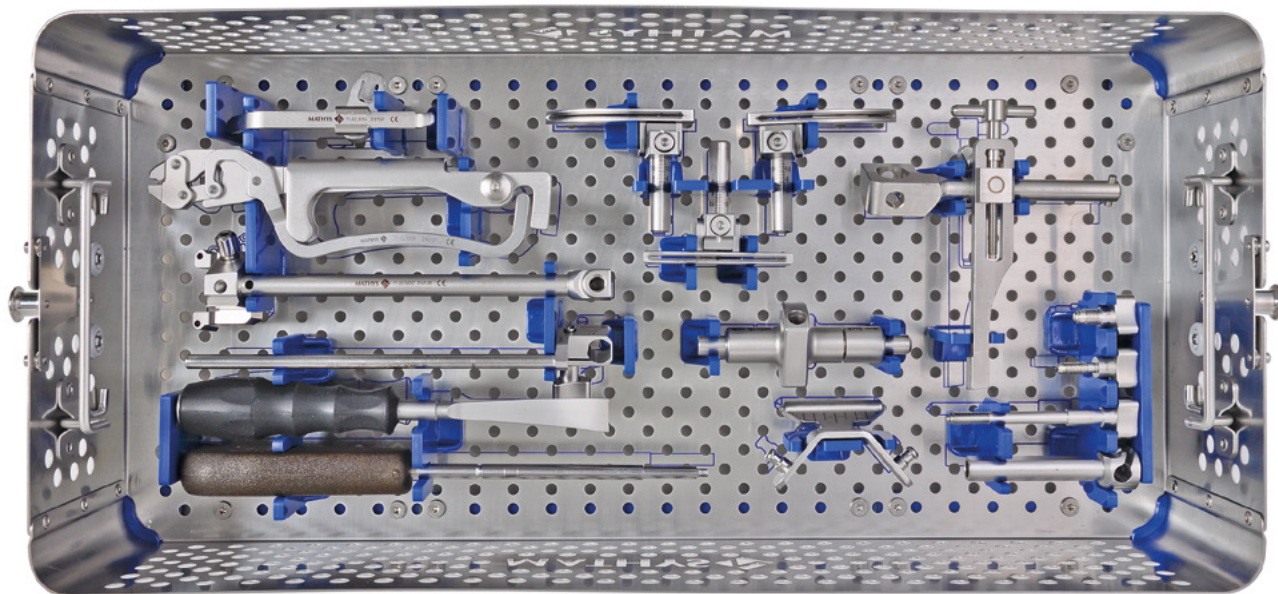

71.34.0546 **balanSys Basic Set Tray**



Item no.	Description	Qty.
71.02.3096	balanSys tibial stylus	1

Item no.	Description	Qty.
71.02.3006	balanSys pliers	1

Item no.	Description	Qty.
71.02.3032	balanSys Trs. aiming device proximal	1

Item no.	Description	Qty.
71.02.3034	balanSys Trs. aiming device distal	1

Item no.	Description	Qty.
71.34.0686	balanSys osteophyte chisel, curved	1

Item no.	Description	Qty.
314.270	Screwdriver, hex., 3.5	1

Item no.	Description	Qty.
71.02.4018	balanSys tibial cutting guide 1.3	1

Item no.	Description	Qty.
71.02.3083	balanSys Trs. tibial cut. guide ri. LIS	1

Item no.	Description	Qty.
71.02.3084	balanSys Trs. tibial cut. guide lft LIS	1

Item no.	Description	Qty.
71.02.3043	balanSys Trs. handle f/intramedullary rod	1

Item no.	Description	Qty.
71.02.3035	balanSys Trs. ankle holder	1



Item no.	Description	Qty.
77.02.0041	balanSys Trs. connecting screw	1


Item no.	Description	Qty.
77.02.0019	balanSys Trs. screw	1

Item no.	Description	Qty.
77.02.0043	balanSys Trs. locking bolt	1

Item no.	Description	Qty.
71.02.3036	balanSys Trs. distal connector	1

Item no.	Description	Qty.
71.02.3041	balanSys Trs. intramedullary shackle	1

balanSys Basic Set 71.34.9000A

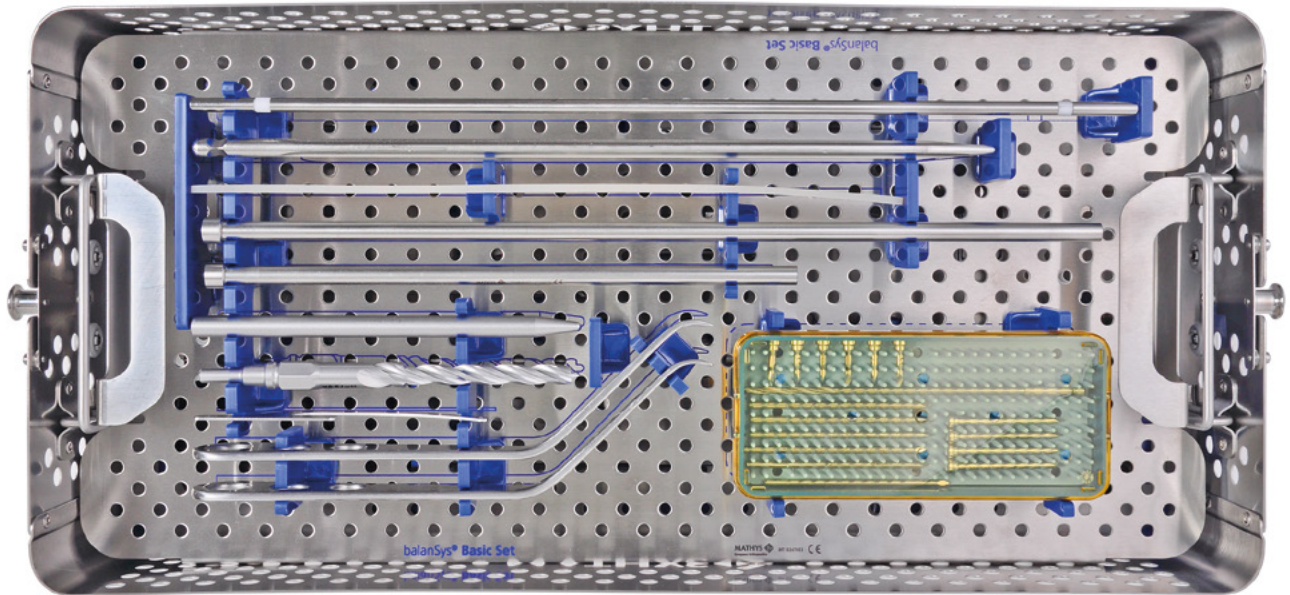

71.34.0547 balanSys Basic Set Tray Insert



Item no.	Description	Qty.
70.04.0109	balanSys alignment rod centre piece	1



Item no.	Description	Qty.
71.02.3042	balanSys Trs. intramedullary rod	1



Item no.	Description	Qty.
71.02.1005	balanSys Trs. rubber band 3x25x300	1



Item no.	Description	Qty.
70.04.0111	balanSys alignment rod long	1



Item no.	Description	Qty.
70.04.0110	balanSys alignment rod short	1



Item no.	Description	Qty.
71.02.3014	balanSys impaction/extraction rod	1



Item no.	Description	Qty.
71.02.3054	balanSys pin 3.2/80	6



Item no.	Description	Qty.
71.34.0108	balanSys Pin with Head 3.2/20	6



Item no.	Description	Qty.
71.02.3003	Pin with head, Ø3.2/6.5x52	4



Item no.	Description	Qty.
71.34.0100	balanSys drill bit 8.5/11 mm	1



Item no.	Description	Qty.
315.310	AO Drill bit 3.2	1



Item no.	Description	Qty.
71.34.0107	balanSys feeler plate 1.3	1



Item no.	Description	Qty.
71.02.3005	balanSys bone retractor	2



Optional instruments


Item no.	Description	Qty.
71.34.0647	Drill Pin 3.2 mm x 89 mm	6



Item no.	Description	Qty.
71.34.0787	Quick Coupling Square 2.25 (Adapter for Drill Pin)	1

balanSys Basic Trial Set Essential 71.34.9005A

No Picture / 71.34.0203 **balanSys Trial Set Essential No. 1 Lid**

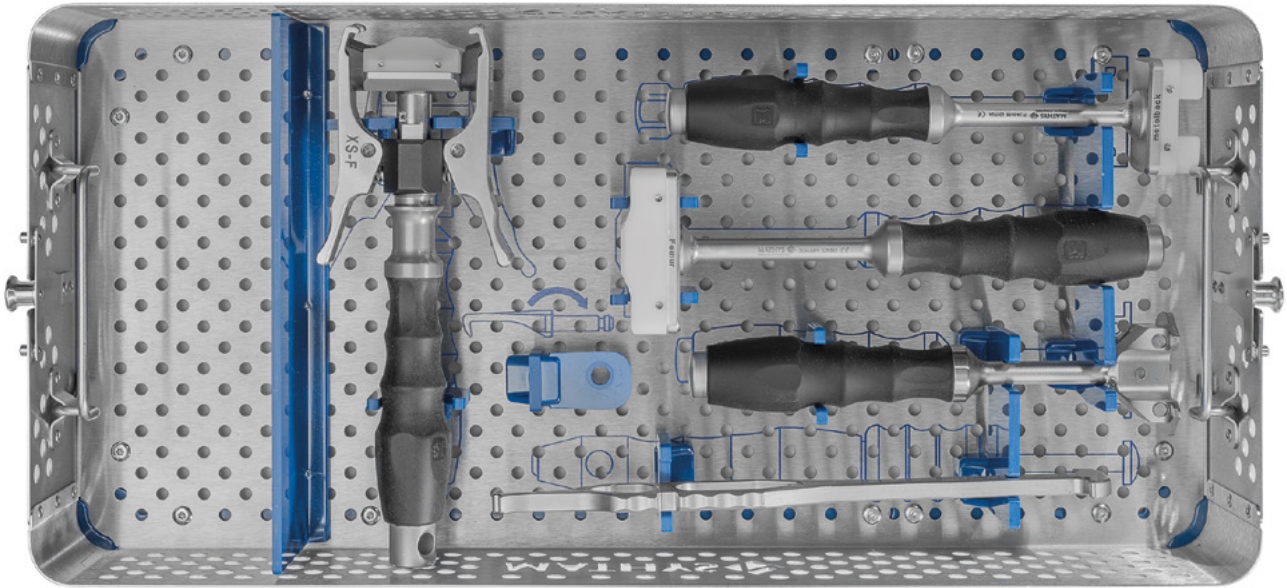

71.34.0204 **balanSys Trial Set Essential No. 1 Tray**



Item no.	Description	Qty.
71.34.0744	balanSys Femur Holder XS-F	1



Item no.	Description	Qty.
71.34.0698	balanSys tibial impactor	1



Item no.	Description	Qty.
71.34.0699	balanSys femoral impactor	1



Item no.	Description	Qty.
71.34.0240	balanSys positioner for tibial plateau	1



Item no.	Description	Qty.
71.34.0788	balanSys Femur Extractor	1




Optional instruments

Item no.	Description	Qty.
71.34.0745	balanSys Femur Holder XS-D	1

balanSys Basic Trial Set Essential 71.34.9005A

No Picture / 71.34.0205 **balanSys Trial Set Essential No.2 Lid**

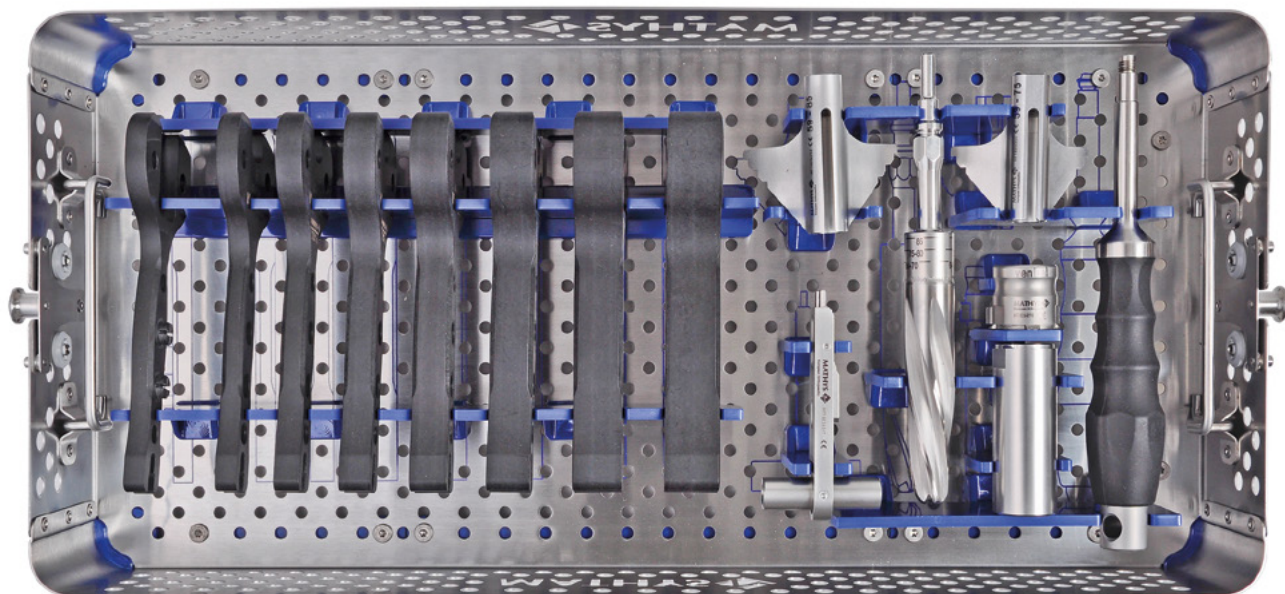

71.34.0206 **balanSys Trial Set Essential No.2 Tray**



Item no.	Description	Qty.
79.02.0651	balanSys Spacer block femur 9	1

Item no.	Description	Qty.
79.02.0640	balanSys Spacer block tibia 8	1
71.34.0947*	balanSys Spacer block tibia 9	1
79.02.0641	balanSys Spacer block tibia 10.5	1
71.34.0948*	balanSys Spacer block tibia 11.5	1
79.02.0642	balanSys Spacer block tibia 13	1
79.02.0643	balanSys Spacer block tibia 15.5	1
79.02.0644	balanSys Spacer block tibia 18	1
79.02.0645	balanSys Spacer block tibia 20.5	1
79.02.0646	balanSys Spacer block tibia 23	1

* balanSys PE Inlays 9mm and 11.5 mm are available in vitamys only.



Item no.	Description	Qty.
71.34.0198	balanSys Fin Chisel 59–85	1
71.34.0199	balanSys Fin Chisel 59–75	1



Item no.	Description	Qty.
71.34.0197	balanSys Chisel Centering Guide	1



Item no.	Description	Qty.
71.34.0196	balanSys holder tibial template	1



Item no.	Description	Qty.
71.34.0200	balanSys Reamer	1




Item no.	Description	Qty.
71.34.0201	balanSys Reamer Guide	1



Item no.	Description	Qty.
71.34.0700	balanSys handle for tib. cutting chisel	1

balanSys LIS 4in1 71.34.9030A

No Picture / 71.34.0162 **balanSys LIS 4in1 Lid**

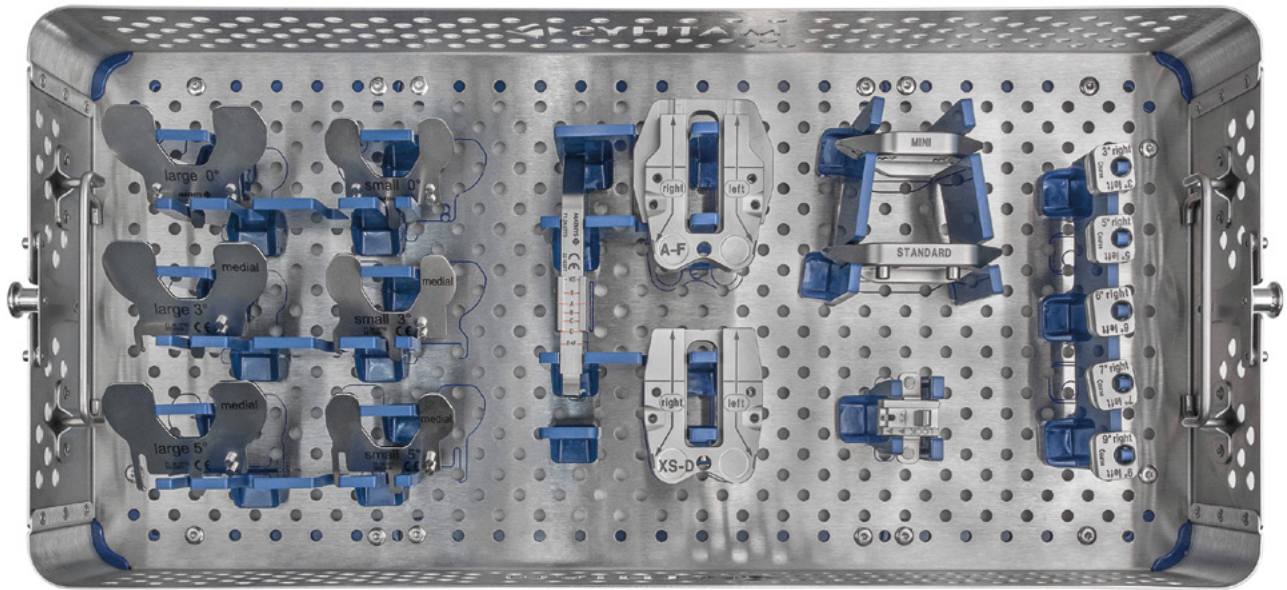

71.34.0163 **balanSys LIS 4in1 Tray**



Item no.	Description	Qty.
71.34.0767	balanSys Fem. Sizing/Rot. Guide LIS XS–D	1

Item no.	Description	Qty.
71.34.0768	balanSys Fem. Sizing/Rot. Guide LIS A–F	1

Item no.	Description	Qty.
71.34.0769	balanSys Femur Rotation Bearing Small 0°	1
71.34.0770	balanSys Femur Rotation Bearing Small 3°	1
71.34.0771	balanSys Femur Rotation Bearing Small 5°	1

Item no.	Description	Qty.
71.34.0772	balanSys Femur Rotation Bearing Large 0°	1
71.34.0773	balanSys Femur Rotation Bearing Large 3°	1
71.34.0774	balanSys Femur Rotation Bearing Large 5°	1

Item no.	Description	Qty.
71.34.0775	balanSys Size Feeler LIS	1

Item no.	Description	Qty.
71.34.0776	balanSys Holder Distal Cutting Guide LIS	1

Item no.	Description	Qty.
71.02.3091	balanSys axis alignment block 3° LIS	1
71.02.3092	balanSys axis alignment block 5° LIS	1
71.34.0777	balanSys axis alignment block 6° LIS	1
71.02.3093	balanSys axis alignment block 7° LIS	1
71.02.3094	balanSys axis alignment block 9° LIS	1

Item no.	Description	Qty.
71.34.0121	balanSys Distal Cutting Guide STANDARD	1


Optional instruments

NOT part of the standard configuration and must be ordered separately:

Item no.	Description	Qty.
71.34.0128	balanSys Distal Cutting Guide MINI	1

balanSys Trial Set CR/UC No. 1

No Picture / 71.34.0241 **balanSys Trial Set CR/UC No. 1 Lid**

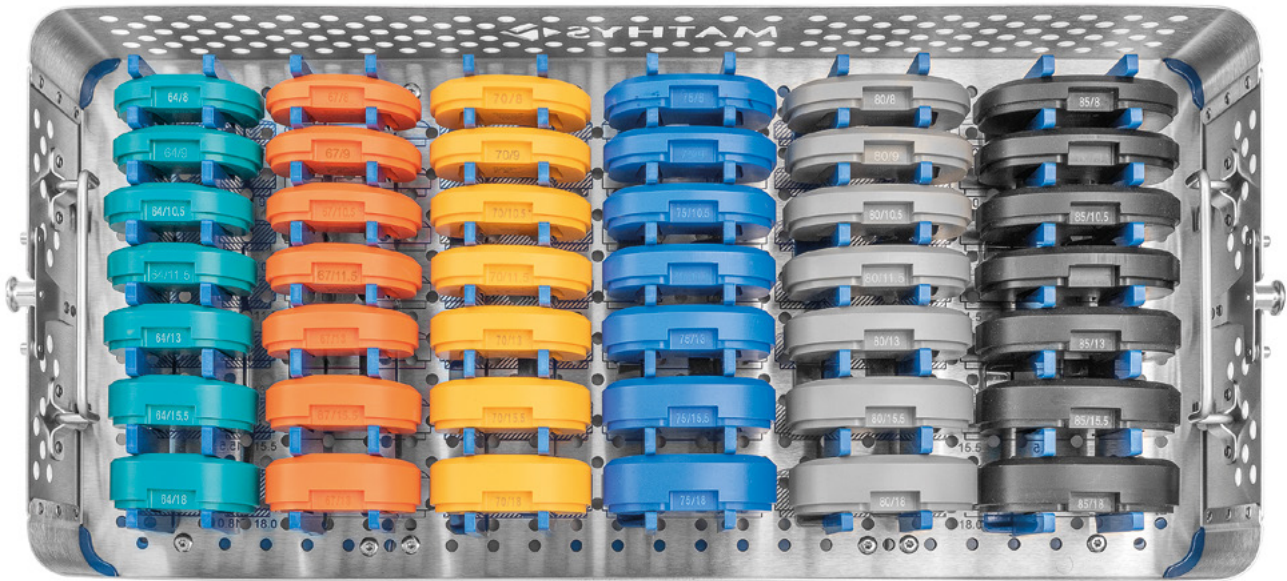

71.34.1075 **balanSys Trial Set 7-CR/UC No. 1 Tray**

balanSys Trial Set CR/UC
Choose your size configuration

Set-no.	Femur	Tibia
71.34.9050A	A/B/C/D/E	64/67/70/75/80/85
71.34.0789A	XS/S/F	59/62



Item no.	Description	Qty.
71.34.0357	balanSys PE Trial Inlay 59/8	1
71.34.0949*	balanSys PE Trial Inlay 59/9	1
71.34.0358	balanSys PE Trial Inlay 59/10.5	1
71.34.0950*	balanSys PE Trial Inlay 59/11.5	1
71.34.0359	balanSys PE Trial Inlay 59/13	1
71.34.0360	balanSys PE Trial Inlay 59/15.5	1
71.34.0361	balanSys PE Trial Inlay 59/18	1
71.34.0362	balanSys PE Trial Inlay 59/20.5	1
71.34.0210	balanSys PE Trial Inlay 62/8	1
71.34.0951*	balanSys PE Trial Inlay 62/9	1
71.34.0211	balanSys PE Trial Inlay 62/10.5	1
71.34.0952*	balanSys PE Trial Inlay 62/11.5	1
71.34.0212	balanSys PE Trial Inlay 62/13	1
71.34.0213	balanSys PE Trial Inlay 62/15.5	1
71.34.0214	balanSys PE Trial Inlay 62/18	1
71.34.0215	balanSys PE Trial Inlay 62/20.5	1
71.34.0216	balanSys PE Trial Inlay 64/8	1
71.34.0953*	balanSys PE Trial Inlay 64/9	1
71.34.0217	balanSys PE Trial Inlay 64/10.5	1
71.34.0954*	balanSys PE Trial Inlay 64/11.5	1
71.34.0218	balanSys PE Trial Inlay 64/13	1
71.34.0219	balanSys PE Trial Inlay 64/15.5	1
71.34.0220	balanSys PE Trial Inlay 64/18	1
71.34.0221	balanSys PE Trial Inlay 64/20.5	1
71.34.0222	balanSys PE Trial Inlay 67/8	1
71.34.0955*	balanSys PE Trial Inlay 67/9	1
71.34.0223	balanSys PE Trial Inlay 67/10.5	1
71.34.0956*	balanSys PE Trial Inlay 67/11.5	1
71.34.0224	balanSys PE Trial Inlay 67/13	1
71.34.0225	balanSys PE Trial Inlay 67/15.5	1
71.34.0226	balanSys PE Trial Inlay 67/18	1
71.34.0227	balanSys PE Trial Inlay 67/20.5	1
71.34.0477	balanSys PE Trial Inlay 70/8	1
71.34.0957*	balanSys PE Trial Inlay 70/9	1
71.34.0478	balanSys PE Trial Inlay 70/10.5	1
71.34.0958*	balanSys PE Trial Inlay 70/11.5	1
71.34.0479	balanSys PE Trial Inlay 70/13	1
71.34.0480	balanSys PE Trial Inlay 70/15.5	1
71.34.0481	balanSys PE Trial Inlay 70/18	1
71.34.0482	balanSys PE Trial Inlay 70/20.5	1


* balanSys PE Inlays 9 mm and 11.5 mm are available in vitamys only.

Item no.	Description	Qty.
71.34.0483	balanSys PE Trial Inlay 75/8	1
71.34.0959*	balanSys PE Trial Inlay 75/9	1
71.34.0484	balanSys PE Trial Inlay 75/10.5	1
71.34.0960*	balanSys PE Trial Inlay 75/11.5	1
71.34.0485	balanSys PE Trial Inlay 75/13	1
71.34.0486	balanSys PE Trial Inlay 75/15.5	1
71.34.0487	balanSys PE Trial Inlay 75/18	1
71.34.0488	balanSys PE Trial Inlay 75/20.5	1
71.34.0228	balanSys PE Trial Inlay 80/8	1
71.34.0961*	balanSys PE Trial Inlay 80/9	1
71.34.0229	balanSys PE Trial Inlay 80/10.5	1
71.34.0962*	balanSys PE Trial Inlay 80/11.5	1
71.34.0230	balanSys PE Trial Inlay 80/13	1
71.34.0231	balanSys PE Trial Inlay 80/15.5	1
71.34.0232	balanSys PE Trial Inlay 80/18	1
71.34.0233	balanSys PE Trial Inlay 80/20.5	1
71.34.0234	balanSys PE Trial Inlay 85/8	1
71.34.0963*	balanSys PE Trial Inlay 85/9	1
71.34.0235	balanSys PE Trial Inlay 85/10.5	1
71.34.0964*	balanSys PE Trial Inlay 85/11.5	1
71.34.0236	balanSys PE Trial Inlay 85/13	1
71.34.0237	balanSys PE Trial Inlay 85/15.5	1
71.34.0238	balanSys PE Trial Inlay 85/18	1
71.34.0239	balanSys PE Trial Inlay 85/20.5	1

* balanSys PE Inlays 9 mm and 11.5 mm are available in vitamys only.

balanSys Trial Set CR/UC No. 2

No Picture / 71.34.0243 balanSys Trial Set CR/UC No. 2 Lid

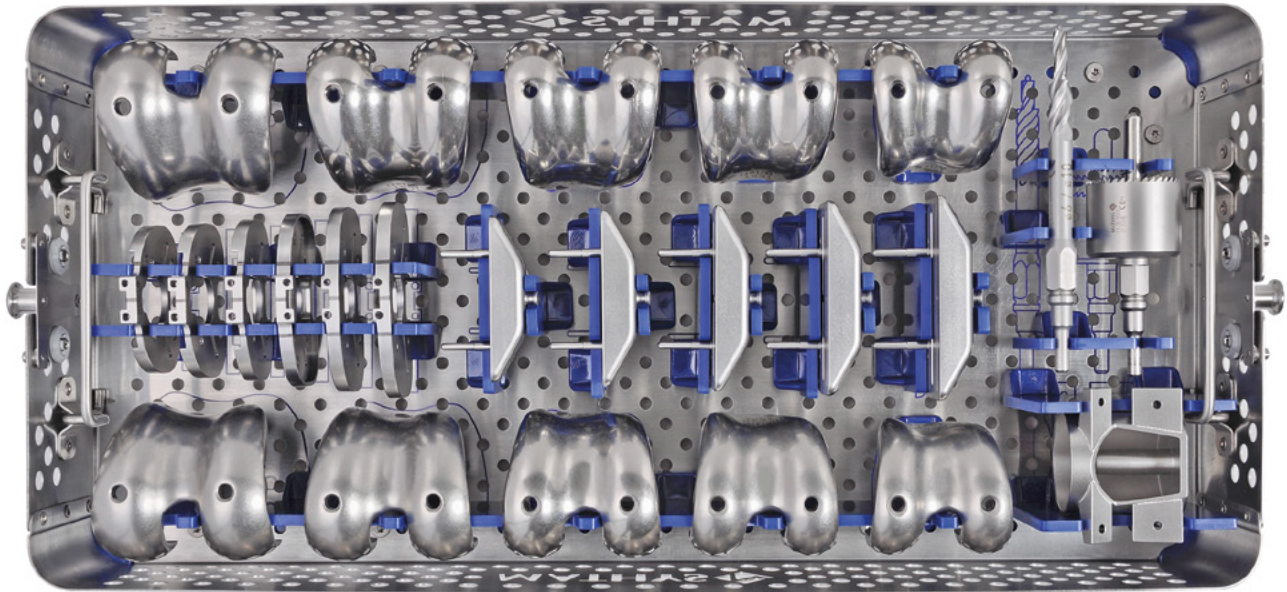

71.34.0244 balanSys Trial Set CR/UC No. 2 Tray



Item no.	Description	Qty.
71.34.0355	balanSys trial femur XS left	1
71.34.0356	balanSys trial femur XS right	1
71.34.0504	balanSys trial femur S left	1
71.34.0505	balanSys trial femur S right	1
71.02.4001	balanSys trial femur A left	1
71.02.4002	balanSys trial femur A right	1
71.02.4301	balanSys trial femur B left	1
71.02.4302	balanSys trial femur B right	1
71.02.4601	balanSys trial femur C left	1
71.02.4602	balanSys trial femur C right	1
71.02.4901	balanSys trial femur D left	1
71.02.4902	balanSys trial femur D right	1
71.02.5201	balanSys trial femur E left	1
71.02.5202	balanSys trial femur E right	1
71.34.0371	balanSys trial femur F left	1
71.34.0372	balanSys trial femur F right	1



Item no.	Description	Qty.
71.34.0353V	balanSys 4in1 Cutting Block STANDARD XS	1
71.34.0122V	balanSys 4in1 Cutting Block STANDARD S	1
71.34.0123V	balanSys 4in1 Cutting Block STANDARD A	1
71.34.0124V	balanSys 4in1 Cutting Block STANDARD B	1
71.34.0125V	balanSys 4in1 Cutting Block STANDARD C	1
71.34.0126V	balanSys 4in1 Cutting Block STANDARD D	1
71.34.0127V	balanSys 4in1 Cutting Block STANDARD E	1
71.34.0370V	balanSys 4in1 Cutting Block STANDARD F	1



Item no.	Description	Qty.
71.34.0535	balanSys CR/PS Tibial Template 59	1
71.34.0536	balanSys CR/PS Tibial Template 62	1
71.34.0537	balanSys CR/PS Tibial Template 64	1
71.34.0538	balanSys CR/PS Tibial Template 67	1
71.34.0539	balanSys CR/PS Tibial Template 70	1
71.34.0540	balanSys CR/PS Tibial Template 75	1
71.34.0541	balanSys CR/PS Tibial Template 80	1
71.34.0542	balanSys CR/PS Tibial Template 85	1



Item no.	Description	Qty.
71.02.3023	balanSys trochlea reamer	1



Item no.	Description	Qty.
71.02.3024	balanSys trochlea bushing	1



Item no.	Description	Qty.
71.34.0023	balanSys Drill Bit with stop 6	1

Optional instruments


NOT part of the standard configuration and must be ordered separately:



Item no.	Description	Qty.
71.34.0363	balanSys 4in1 Cutting Block MINI XS	1
71.34.0129	balanSys 4in1 Cutting Block MINI S	1
71.34.0130	balanSys 4in1 Cutting Block MINI A	1
71.34.0131	balanSys 4in1 Cutting Block MINI B	1
71.34.0132	balanSys 4in1 Cutting Block MINI C	1
71.34.0133	balanSys 4in1 Cutting Block MINI D	1
71.34.0134	balanSys 4in1 Cutting Block MINI E	1
71.34.0373	balanSys 4in1 Cutting Block MINI F	1

balanSys Trial Set PS

No Picture / 71.34.0286 **balanSys trial set PS no. 1 Lid**

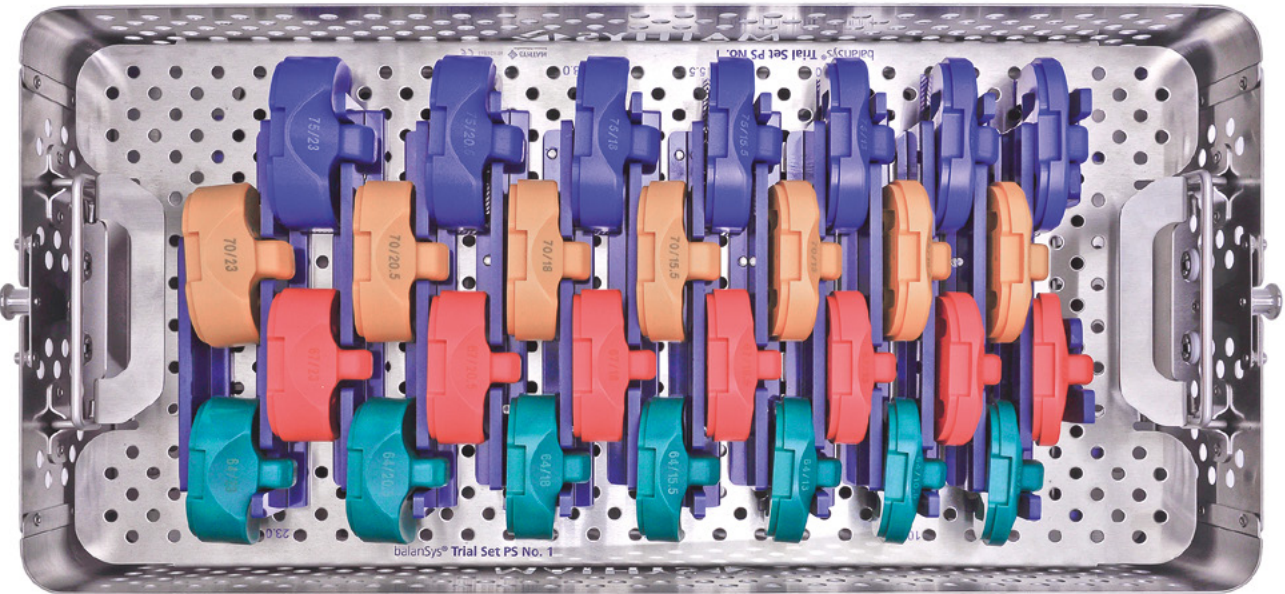

71.34.0288 **balanSys Trial Set PS No. 1 Tray Insert**

balanSys Trial Set PS
Choose your size configuration

Set-no.	Femur	Tibia
71.34.9070A	A/B/C/D/E	64/67/70/75/80/85
71.34.0790A	XS/S/F	59/62




Item no.	Description	Qty.
71.34.0384	balanSys PS PE Trial Inlay 59/8	1
71.34.0965*	balanSys PS PE Trial Inlay 59/9	1
71.34.0385	balanSys PS PE Trial Inlay 59/10.5	1
71.34.0966*	balanSys PS PE Trial Inlay 59/11.5	1
71.34.0386	balanSys PS PE Trial Inlay 59/13	1
71.34.0387	balanSys PS PE Trial Inlay 59/15.5	1
71.34.0388	balanSys PS PE Trial Inlay 59/18	1
71.34.0389	balanSys PS PE Trial Inlay 59/20.5	1
71.34.0503	balanSys PS PE Trial Inlay 59/23	1
71.34.0249	balanSys PS PE Trial Inlay 62/8	1
71.34.0967*	balanSys PS PE Trial Inlay 62/9	1
71.34.0250	balanSys PS PE Trial Inlay 62/10.5	1
71.34.0968*	balanSys PS PE Trial Inlay 62/11.5	1
71.34.0251	balanSys PS PE Trial Inlay 62/13	1
71.34.0252	balanSys PS PE Trial Inlay 62/15.5	1
71.34.0253	balanSys PS PE Trial Inlay 62/18	1
71.34.0254	balanSys PS PE Trial Inlay 62/20.5	1
71.34.0255	balanSys PS PE Trial Inlay 62/23	1
71.34.0256	balanSys PS PE Trial Inlay 64/8	1
71.34.0969*	balanSys PS PE Trial Inlay 64/9	1
71.34.0257	balanSys PS PE Trial Inlay 64/10.5	1
71.34.0970*	balanSys PS PE Trial Inlay 64/11.5	1
71.34.0258	balanSys PS PE Trial Inlay 64/13	1
71.34.0259	balanSys PS PE Trial Inlay 64/15.5	1
71.34.0260	balanSys PS PE Trial Inlay 64/18	1
71.34.0261	balanSys PS PE Trial Inlay 64/20.5	1
71.34.0262	balanSys PS PE Trial Inlay 64/23	1
71.34.0263	balanSys PS PE Trial Inlay 67/8	1
71.34.0971*	balanSys PS PE Trial Inlay 67/9	1
71.34.0264	balanSys PS PE Trial Inlay 67/10.5	1
71.34.0972*	balanSys PS PE Trial Inlay 67/11.5	1
71.34.0265	balanSys PS PE Trial Inlay 67/13	1
71.34.0266	balanSys PS PE Trial Inlay 67/15.5	1
71.34.0267	balanSys PS PE Trial Inlay 67/18	1
71.34.0268	balanSys PS PE Trial Inlay 67/20.5	1
71.34.0269	balanSys PS PE Trial Inlay 67/23	1

* balanSys PE Inlays 9 mm and 11.5 mm are available in vitamys only.

Item no.	Description	Qty.
71.34.0489	balanSys PS PE Trial Inlay 70/8	1
71.34.0973*	balanSys PS PE Trial Inlay 70/9	1
71.34.0490	balanSys PS PE Trial Inlay 70/10.5	1
71.34.0974*	balanSys PS PE Trial Inlay 70/11.5	1
71.34.0491	balanSys PS PE Trial Inlay 70/13	1
71.34.0492	balanSys PS PE Trial Inlay 70/15.5	1
71.34.0493	balanSys PS PE Trial Inlay 70/18	1
71.34.0494	balanSys PS PE Trial Inlay 70/20.5	1
71.34.0495	balanSys PS PE Trial Inlay 70/23	1
71.34.0496	balanSys PS PE Trial Inlay 75/8	1
71.34.0975*	balanSys PS PE Trial Inlay 75/9	1
71.34.0497	balanSys PS PE Trial Inlay 75/10.5	1
71.34.0976*	balanSys PS PE Trial Inlay 75/11.5	1
71.34.0498	balanSys PS PE Trial Inlay 75/13	1
71.34.0499	balanSys PS PE Trial Inlay 75/15.5	1
71.34.0500	balanSys PS PE Trial Inlay 75/18	1
71.34.0501	balanSys PS PE Trial Inlay 75/20.5	1
71.34.0502	balanSys PS PE Trial Inlay 75/23	1
71.34.0270	balanSys PS PE Trial Inlay 80/8	1
71.34.0977*	balanSys PS PE Trial Inlay 80/9	1
71.34.0271	balanSys PS PE Trial Inlay 80/10.5	1
71.34.0978*	balanSys PS PE Trial Inlay 80/11.5	1
71.34.0272	balanSys PS PE Trial Inlay 80/13	1
71.34.0273	balanSys PS PE Trial Inlay 80/15.5	1
71.34.0274	balanSys PS PE Trial Inlay 80/18	1
71.34.0275	balanSys PS PE Trial Inlay 80/20.5	1
71.34.0276	balanSys PS PE Trial Inlay 80/23	1
71.34.0277	balanSys PS PE Trial Inlay 85/8	1
71.34.0979*	balanSys PS PE Trial Inlay 85/9	1
71.34.0278	balanSys PS PE Trial Inlay 85/10.5	1
71.34.0980*	balanSys PS PE Trial Inlay 85/11.5	1
71.34.0279	balanSys PS PE Trial Inlay 85/13	1
71.34.0280	balanSys PS PE Trial Inlay 85/15.5	1
71.34.0281	balanSys PS PE Trial Inlay 85/18	1
71.34.0282	balanSys PS PE Trial Inlay 85/20.5	1
71.34.0283	balanSys PS PE Trial Inlay 85/23	1

* balanSys PE Inlays 9 mm and 11.5 mm are available in vitamys only.

balanSys Trial Set PS

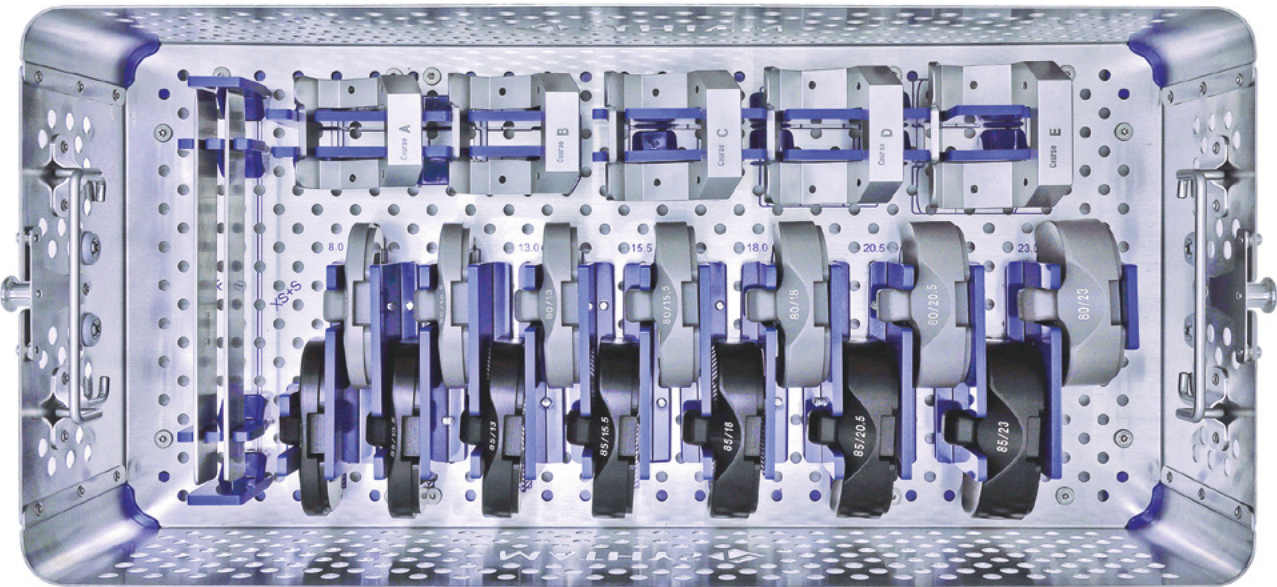

71.34.0287 balanSys Trial Set PS No. 1 Tray



Item no.	Description	Qty.
71.34.0690	balanSys Chisel 22 mm XS/S	1




Item no.	Description	Qty.
71.34.0691	balanSys Chisel 25 mm A–F	1



Item no.	Description	Qty.
71.34.0390	balanSys PS femur box cutting guide XS	1
71.34.0284	balanSys PS femur box cutting guide S	1
79.02.0265	balanSys PS femur box cutting guide A	1
79.02.0266	balanSys PS femur box cutting guide B	1
79.02.0267	balanSys PS femur box cutting guide C	1
79.02.0268	balanSys PS femur box cutting guide D	1
79.02.0269	balanSys PS femur box cutting guide E	1
71.34.0401	balanSys PS femur box cutting guide F	1

balanSys Trial Set PS

No Picture / 71.34.0289 **balanSys Trial Set PS No. 2 Lid**

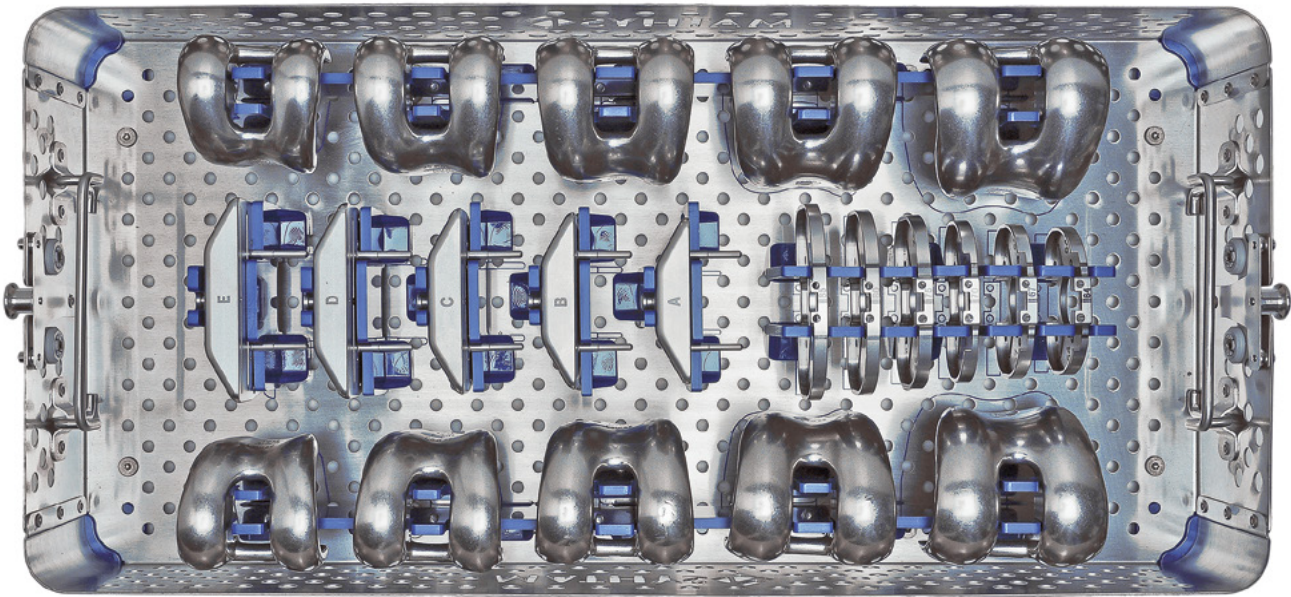

71.34.0290 **balanSys Trial Set PS No. 2 Tray**



Item no.	Description	Qty.
71.34.0383	balanSys PS trial femur XS right	1
71.34.0382	balanSys PS trial femur XS left	1
71.34.0248	balanSys PS trial femur S right	1
71.34.0247	balanSys PS trial femur S left	1
79.02.0040	balanSys PS trial femur A right	1
79.02.0041	balanSys PS trial femur A left	1
79.02.0042	balanSys PS trial femur B right	1
79.02.0043	balanSys PS trial femur B left	1
79.02.0044	balanSys PS trial femur C right	1
79.02.0045	balanSys PS trial femur C left	1
79.02.0046	balanSys PS trial femur D right	1
79.02.0047	balanSys PS trial femur D left	1
79.02.0048	balanSys PS trial femur E right	1
79.02.0049	balanSys PS trial femur E left	1
71.34.0400	balanSys PS trial femur F right	1
71.34.0399	balanSys PS trial femur F left	1



Item no.	Description	Qty.
71.34.0353V	balanSys 4in1 Cutting Block STANDARD XS	1
71.34.0122V	balanSys 4in1 Cutting Block STANDARD S	1
71.34.0123V	balanSys 4in1 Cutting Block STANDARD A	1
71.34.0124V	balanSys 4in1 Cutting Block STANDARD B	1
71.34.0125V	balanSys 4in1 Cutting Block STANDARD C	1
71.34.0126V	balanSys 4in1 Cutting Block STANDARD D	1
71.34.0127V	balanSys 4in1 Cutting Block STANDARD E	1
71.34.0370V	balanSys 4in1 Cutting Block STANDARD F	1



Item no.	Description	Qty.
71.34.0535	balanSys CR/PS Tibial Template 59	1
71.34.0536	balanSys CR/PS Tibial Template 62	1
71.34.0537	balanSys CR/PS Tibial Template 64	1
71.34.0538	balanSys CR/PS Tibial Template 67	1
71.34.0539	balanSys CR/PS Tibial Template 70	1
71.34.0540	balanSys CR/PS Tibial Template 75	1
71.34.0541	balanSys CR/PS Tibial Template 80	1
71.34.0542	balanSys CR/PS Tibial Template 85	1

Optional instruments


NOT part of the standard configuration and must be ordered separately:



Item no.	Description	Qty.
71.34.0363	balanSys 4in1 Cutting Block MINI XS	1
71.34.0129	balanSys 4in1 Cutting Block MINI S	1
71.34.0130	balanSys 4in1 Cutting Block MINI A	1
71.34.0131	balanSys 4in1 Cutting Block MINI B	1
71.34.0132	balanSys 4in1 Cutting Block MINI C	1
71.34.0133	balanSys 4in1 Cutting Block MINI D	1
71.34.0134	balanSys 4in1 Cutting Block MINI E	1
71.34.0373	balanSys 4in1 Cutting Block MINI F	1

balanSys Trial Set RP 71.34.9060A (optional)

No Picture / 71.34.1056 **leggera Set Lid**

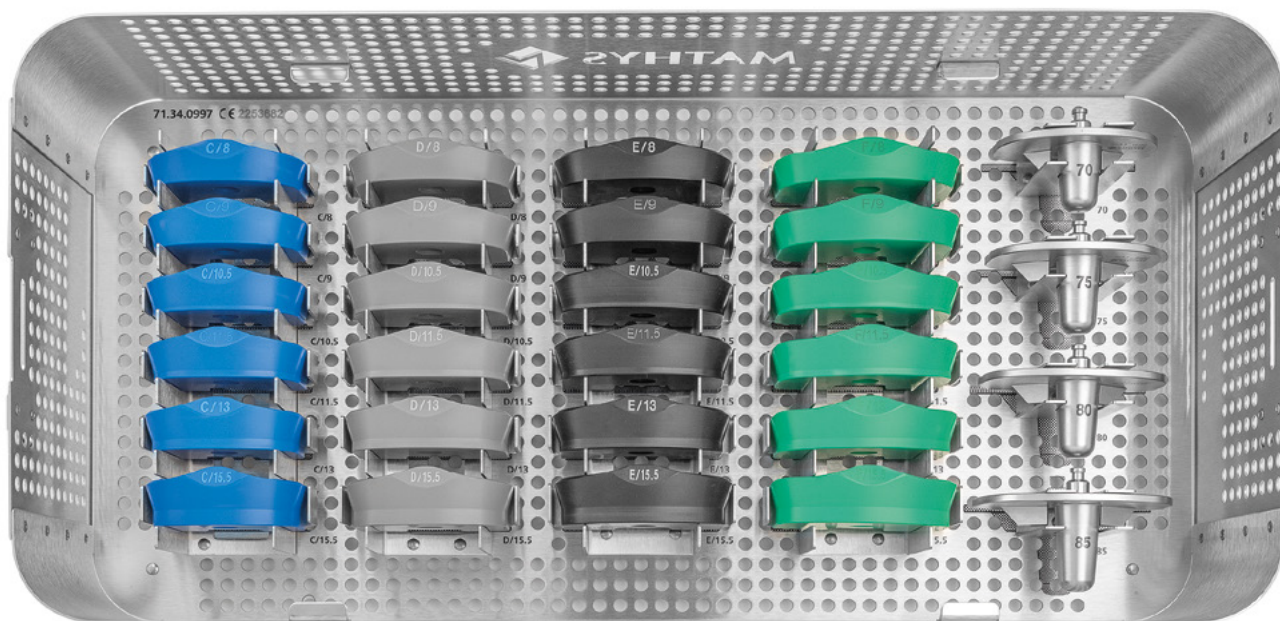

71.34.0997 **balanSys Trial Set 6-RP Tray**




Item no.	Description	Qty.
71.34.0297	balanSys RP Trial Tibial Plateau 70	1
71.34.0298	balanSys RP Trial Tibial Plateau 75	1
71.34.0299	balanSys RP Trial Tibial Plateau 80	1
71.34.0300	balanSys RP Trial Tibial Plateau 85	1

Item no.	Description	Qty.
71.34.0574	balanSys RP PE Trial Inlay C/8	1
71.34.0989*	balanSys RP PE Trial Inlay C/9	1
71.34.0575	balanSys RP PE Trial Inlay C/10.5	1
71.34.0990*	balanSys RP PE Trial Inlay C/11.5	1
71.34.0576	balanSys RP PE Trial Inlay C/13	1
71.34.0577	balanSys RP PE Trial Inlay C/15.5	1
71.34.0580	balanSys RP PE Trial Inlay D/8	1
71.34.0991*	balanSys RP PE Trial Inlay D/9	1
71.34.0581	balanSys RP PE Trial Inlay D/10.5	1
71.34.0992*	balanSys RP PE Trial Inlay D/11.5	1
71.34.0582	balanSys RP PE Trial Inlay D/13	1
71.34.0583	balanSys RP PE Trial Inlay D/15.5	1
71.34.0586	balanSys RP PE Trial Inlay E/8	1
71.34.0993*	balanSys RP PE Trial Inlay E/9	1
71.34.0587	balanSys RP PE Trial Inlay E/10.5	1
71.34.0994*	balanSys RP PE Trial Inlay E/11.5	1
71.34.0588	balanSys RP PE Trial Inlay E/13	1
71.34.0589	balanSys RP PE Trial Inlay E/15.5	1
71.34.0429	balanSys RP PE Trial Inlay F/8	1
71.34.0995*	balanSys RP PE Trial Inlay F/9	1
71.34.0430	balanSys RP PE Trial Inlay F/10.5	1
71.34.0996*	balanSys RP PE Trial Inlay F/11.5	1
71.34.0431	balanSys RP PE Trial Inlay F/13	1
71.34.0432	balanSys RP PE Trial Inlay F/15.5	1

* balanSys PE Inlays 9 mm and 11.5 mm are available in vitamys only.

balanSys Trial Set RP 71.34.9060A (optional)

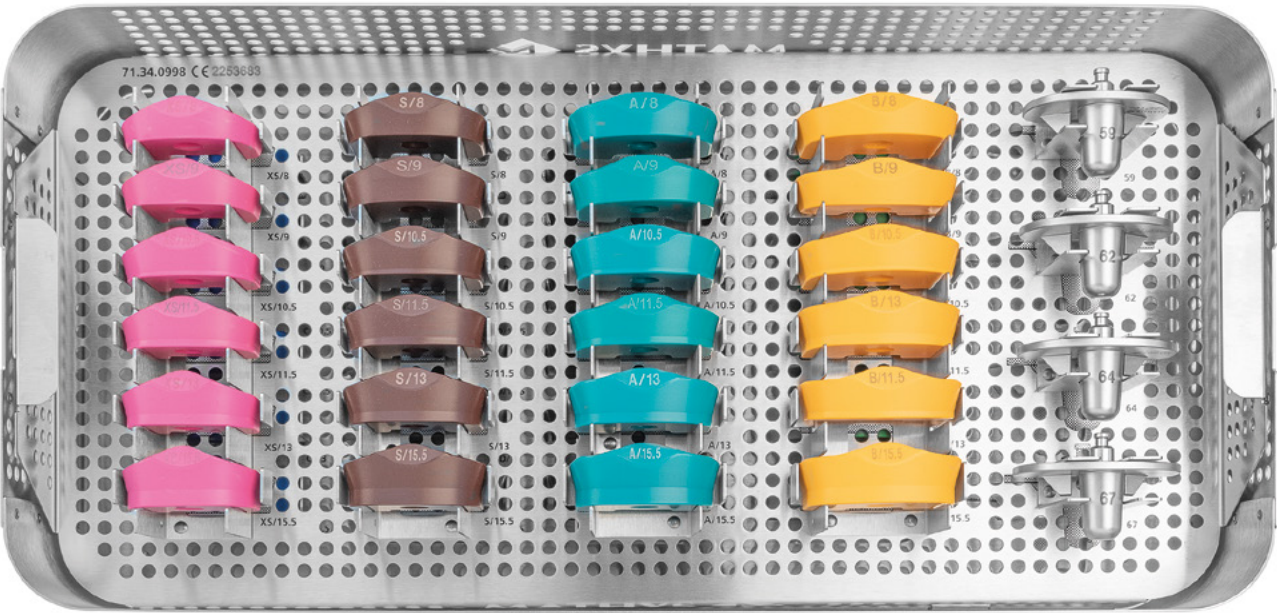

71.34.0998 **balanSys Trial Set 6-RP Insert**



Item no.	Description	Qty.
71.34.0418	balanSys RP Trial Tibial Plateau 59	1
71.34.0294	balanSys RP Trial Tibial Plateau 62	1
71.34.0295	balanSys RP Trial Tibial Plateau 64	1
71.34.0296	balanSys RP Trial Tibial Plateau 67	1

Item no.	Description	Qty.
71.34.0413	balanSys RP PE Trial Inlay XS/8	1
71.34.0981*	balanSys RP PE Trial Inlay XS/9	1
71.34.0414	balanSys RP PE Trial Inlay XS/10.5	1
71.34.0982*	balanSys RP PE Trial Inlay XS/11.5	1
71.34.0415	balanSys RP PE Trial Inlay XS/13	1
71.34.0416	balanSys RP PE Trial Inlay XS/15.5	1
71.34.0301	balanSys RP PE Trial Inlay S/8	1
71.34.0983*	balanSys RP PE Trial Inlay S/9	1
71.34.0302	balanSys RP PE Trial Inlay S/10.5	1
71.34.0984*	balanSys RP PE Trial Inlay S/11.5	1
71.34.0303	balanSys RP PE Trial Inlay S/13	1
71.34.0304	balanSys RP PE Trial Inlay S/15.5	1
71.34.0562	balanSys RP PE Trial Inlay A/8	1
71.34.0985*	balanSys RP PE Trial Inlay A/9	1
71.34.0563	balanSys RP PE Trial Inlay A/10.5	1
71.34.0986*	balanSys RP PE Trial Inlay A/11.5	1
71.34.0564	balanSys RP PE Trial Inlay A/13	1
71.34.0565	balanSys RP PE Trial Inlay A/15.5	1
71.34.0568	balanSys RP PE Trial Inlay B/8	1
71.34.0987*	balanSys RP PE Trial Inlay B/9	1
71.34.0569	balanSys RP PE Trial Inlay B/10.5	1
71.34.0988*	balanSys RP PE Trial Inlay B/11.5	1
71.34.0570	balanSys RP PE Trial Inlay B/13	1
71.34.0571	balanSys RP PE Trial Inlay B/15.5	1

* balanSys PE Inlays 9 mm and 11.5 mm are available in vitamys only.

balanSys 3-Peg Patella FLAT 71.34.0080A


No Picture / 71.34.0082

balanSys Lid 3-Peg Patella FLAT

balanSys 3-Peg Patella Standard 71.34.0081A


No Picture / 71.34.0084

balanSys Lid 3-Peg Patella STANDARD

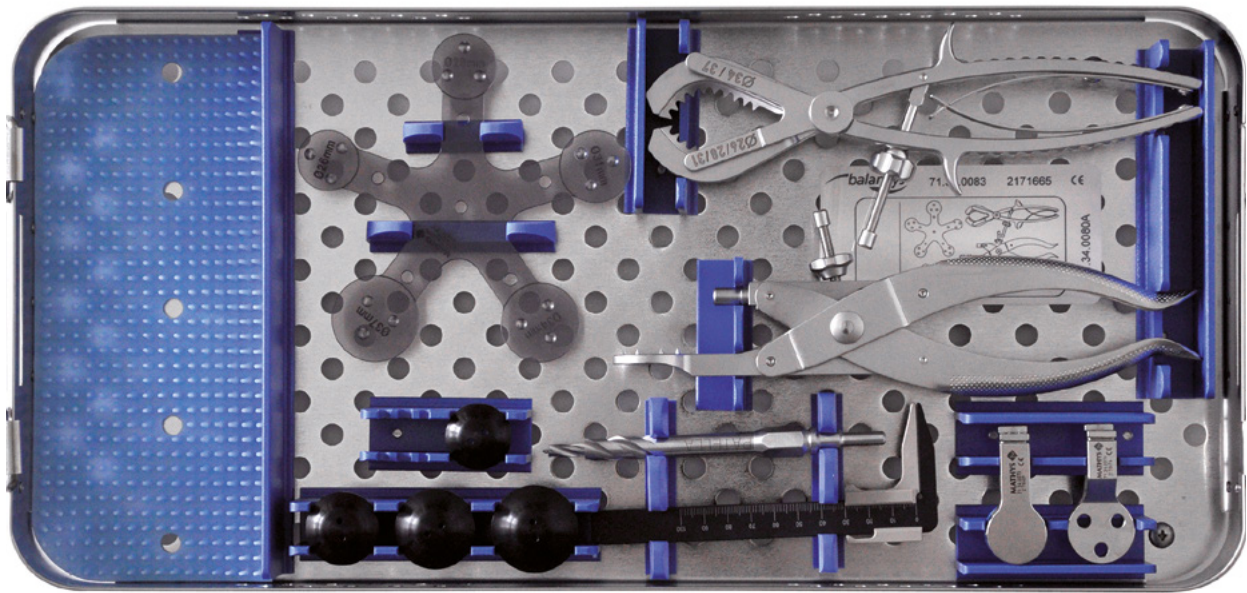


71.34.0083 balanSys Tray f / patella 3 pegs FLAT

71.34.0085 balanSys Tray f / patella 3 pegs STANDARD

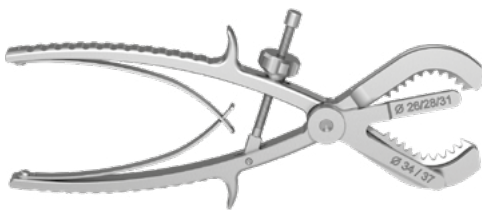

Item no.	Description
71.34.0071	balanSys Patella resection pliers flat

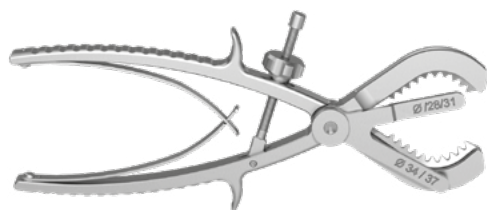

Item no.	Description
71.34.0070	balanSys Patella resection pliers raised




Item no.	Description
71.34.0708	balanSys Trial patella 3 pegs flat 26
71.34.0075	balanSys Trial patella 3 pegs flat 28
71.34.0076	balanSys Trial patella 3 pegs flat 31
71.34.0077	balanSys Trial patella 3 pegs flat 34
71.34.0078	balanSys Trial patella 3 pegs flat 37



Item no.	Description
71.02.3063	balanSys trial patella 3-peg 28
71.02.3064	balanSys trial patella 3-peg 31
71.02.3065	balanSys trial patella 3-peg 34
71.02.3066	balanSys trial patella 3-peg 37



Item no.	Description	Qty.
71.02.2201	balanSys Patella universal pliers	1

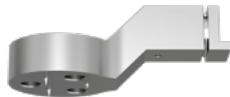


Item no.	Description	Qty.
71.34.0074	balanSys Patella drill guide to pliers	1



Item no.	Description	Qty.
71.34.0073	balanSys Patella cementing aid to pliers	1




Item no.	Description	Qty.
71.02.3061	Drill bit 5.5	1

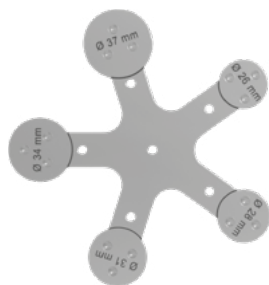

Optional instruments

NOT part of the standard configuration and must be ordered separately:

Item no.	Description	Qty.
71.34.0079	balanSys patella sizing guide	1

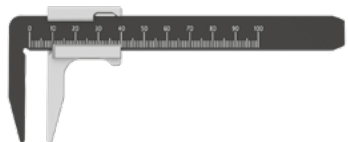


Item no.	Description	Qty.
71.02.3002	balanSys Patella calliper	1

Appendix


6 – Item numbers of the balanSys measuring templates

balanSys BICON Knee System 330.030.034

Suitable for CR, UC and RP

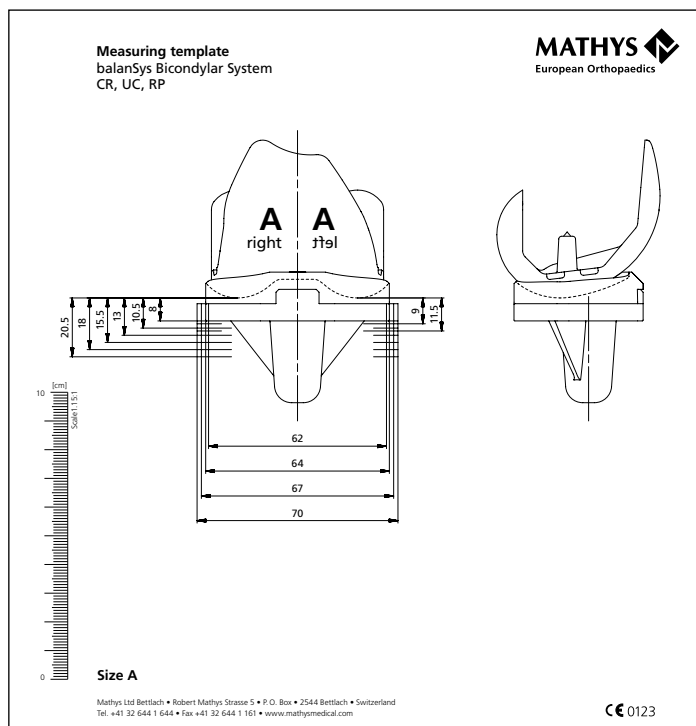

balanSys PS Knee System 330.030.035

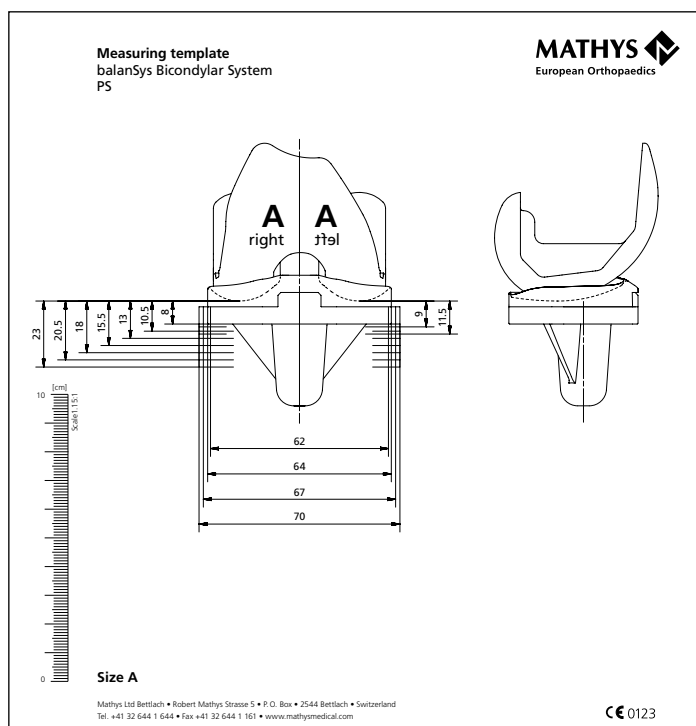

Symbols and abbreviations



Manufacturer



Correct



Incorrect



Caution

CR Cruciate Retaining

UC Ultra Congruent

PS Posterior Stabilized

RP Rotating Platform

ACL Anterior Cruciate Ligament

PCL Posterior Cruciate Ligament

MCL Medial Collateral Ligament

LCL Lateral Collateral Ligament

TRS Tibia Reference System

IFU Instruction For Use

Notes

[illegible]

Australia	Mathys Orthopaedics Pty Ltd Artarmon, NSW 2064 Tel: +61 2 9417 9200 info.au@mathysmedical.com	Italy	Mathys Ortopedia S.r.l. 20141 Milan Tel: +39 02 4959 8085 info.it@mathysmedical.com
Austria	Mathys Orthopädie GmbH 2351 Wiener Neudorf Tel: +43 2236 860 999 info.at@mathysmedical.com	Japan	Mathys KK Tokyo 108-0075 Tel: +81 3 3474 6900 info.jp@mathysmedical.com
Belgium	Mathys Orthopaedics Belux N.V.-S.A. 3001 Leuven Tel: +32 16 38 81 20 info.be@mathysmedical.com	New Zealand	Mathys Ltd. Auckland Tel: +64 9 478 39 00 info.nz@mathysmedical.com
France	Mathys Orthopédie S.A.S 63360 Gerzat Tel: +33 4 73 23 95 95 info.fr@mathysmedical.com	Netherlands	Mathys Orthopaedics B.V. 3001 Leuven Tel: +31 88 1300 500 info.nl@mathysmedical.com
Germany	Mathys Orthopädie GmbH «Centre of Excellence Sales» Bochum 44809 Bochum Tel: +49 234 588 59 0 sales.de@mathysmedical.com «Centre of Excellence Ceramics» Mörsdorf 07646 Mörsdorf/Thür. Tel: +49 364 284 94 0 info.de@mathysmedical.com «Centre of Excellence Production» Hermsdorf 07629 Hermsdorf Tel: +49 364 284 94 110 info.de@mathysmedical.com	P. R. China	Mathys (Shanghai) Medical Device Trading Co., Ltd Shanghai, 200041 Tel: +86 21 6170 2655 info.cn@mathysmedical.com
		Switzerland	Mathys (Schweiz) GmbH 2544 Bettlach Tel: +41 32 644 1 458 info@mathysmedical.com
		United Kingdom	Mathys Orthopaedics Ltd Alton, Hampshire GU34 2QL Tel: +44 8450 580 938 info.uk@mathysmedical.com

Local Marketing Partners in over 30 countries worldwide ...

