


*Distributed by*

**MATHYS**   
European Orthopaedics

Sterile Sawblades  
**Product information**

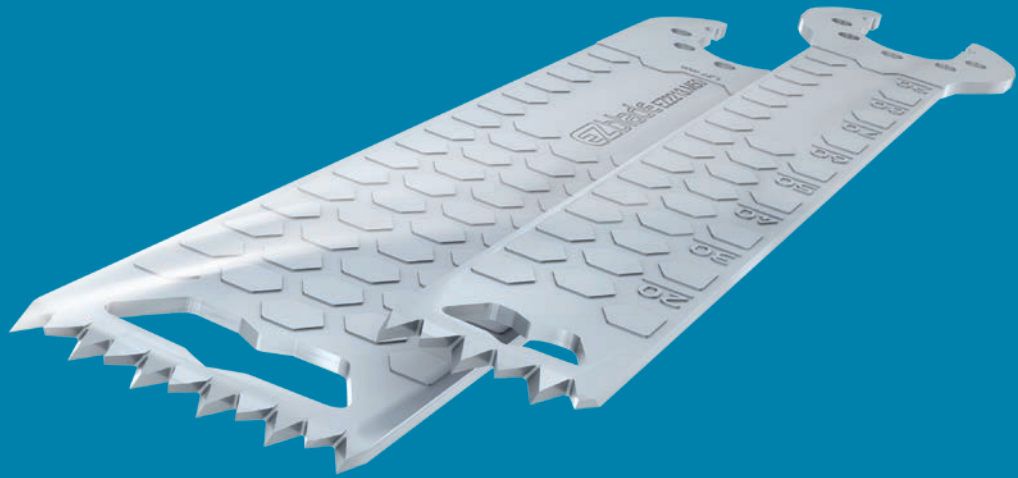

*Manufacturer*

**Gebr. Brasseler  
GmbH & Co. KG**


# EZ Blades – Premium Quality



## The Blades

The enhanced blade geometry with specially ground clearance allows even better evacuation of bone chips and consequently easier penetration of the bone. The result:


- Less heat generation
- Reduced risk of bone necrosis
- A longer service life

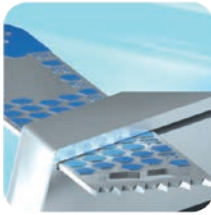

## The Surface

Provided with a honeycomb structure, the ingenious surface structure

- Decreases friction in the saw guide
- Reduces material abrasion

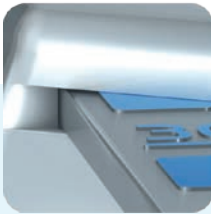

## The Edges

The durability of a saw blade depends not only on the sharpness of its blade, but also on the extent of burr formation at its outer edges. The formation of burr during the sawing process cannot be completely avoided, but in the EZ Blades, this is greatly reduced thanks to their bevelled outer edges.

- More comfort during sawing in the cutting blocks

## Connections available



### Connection

### Compatible Systems

#### **DePuy Synthes**

Compact Air Drive II • Power Drive  
Air Power Line • Batterie Power Line, TRS

#### **Zimmer**

Zimmer Universal Power System



#### **Stryker**

System EHD • System 2000 • System 4 •  
System 5 • System 6 • System 7



#### **Aesculap**

Acculan/Comfort  
(not compatible with Acculan 3Ti)




#### **ConMed Hall**

M Power • M Power 2 • PowerPro Max • PowerPro  
+ PowerPro Electric II • Power Pro Pneumatic •  
Hall Serie 4 • Hall Serie 3 • Versipower Plus • Hall 50

# Sawblades for balanSys® Knee System


The following sawblades are compatible with the **balanSys** instruments:

## Premium Sawblades

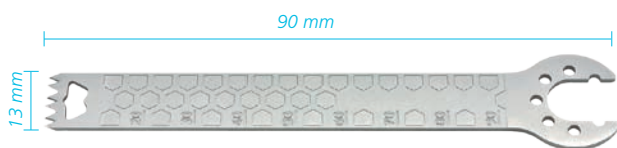

**Sawblade EZ sterile 90x13x1.27**

Item no.	Connection	Dimension
71.34.0622	DePuy Synthes	90x13x1.27

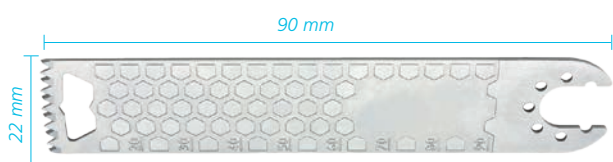

**Sawblade EZ sterile 90x22x1.27**

Item no.	Connection	Dimension
71.34.0618	DePuy Synthes	90x22x1.27
71.34.0693	ConMed Hall	90x22x1.27
71.34.0778	Stryker	90x22x1.27

## Standard Sawblades (Single use)



**Sawblade sterile 90x22x1.27**

Item no.	Connection	Dimension
71.02.3101	DePuy Synthes	90x22x1.27
71.34.0089	ConMed Hall	90x22x1.27
71.34.0090	Aesculap	90x22x1.27
71.34.0753	Stryker	90x22x1.27



**Sawblade sterile 90 x 13 x 1.27**

Item no.	Connection	Dimension
71.02.3106	DePuy Synthes	90x13x1.27
71.34.0092	ConMed Hall	90x13x1.27
71.34.0093	Aesculap	90x13x1.27
71.34.0752	Stryker	90x13x1.27



**Sawblade sterile MIS 100x17x1.27**

Item no.	Connection	Dimension
71.34.0095	DePuy Synthes	100x17x1.27
71.34.0751	Stryker	100x17x1.27

### Sterile-packed

Our carefully chosen selection of knee sawblades are all sterile and individually packed.

Manufacturer

Gebr. Brasseler  
GmbH & Co. KG

Distributed by

**MATHYS**   
European Orthopaedics

## Sawblades for Affinis<sup>®</sup> Shoulder System


The following sawblades are compatible with the **Affinis** instruments:

### Standard Sawblades (Single use)



**Sawblade sterile 90x22x0.89**

Item no.	Connection	Dimension
71.02.03111	DePuy Synthes	90x22x0.89

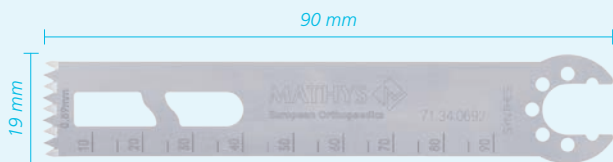

**Sawblade sterile 90x19x0.89**

Item no.	Connection	Dimension
71.34.0692	DePuy Synthes	90x19x0.89

The shoulder sawblades are all sterile and individually packed.