


Manufacturer

Enztec Limited

Distributed by

MATHYS 
European Orthopaedics

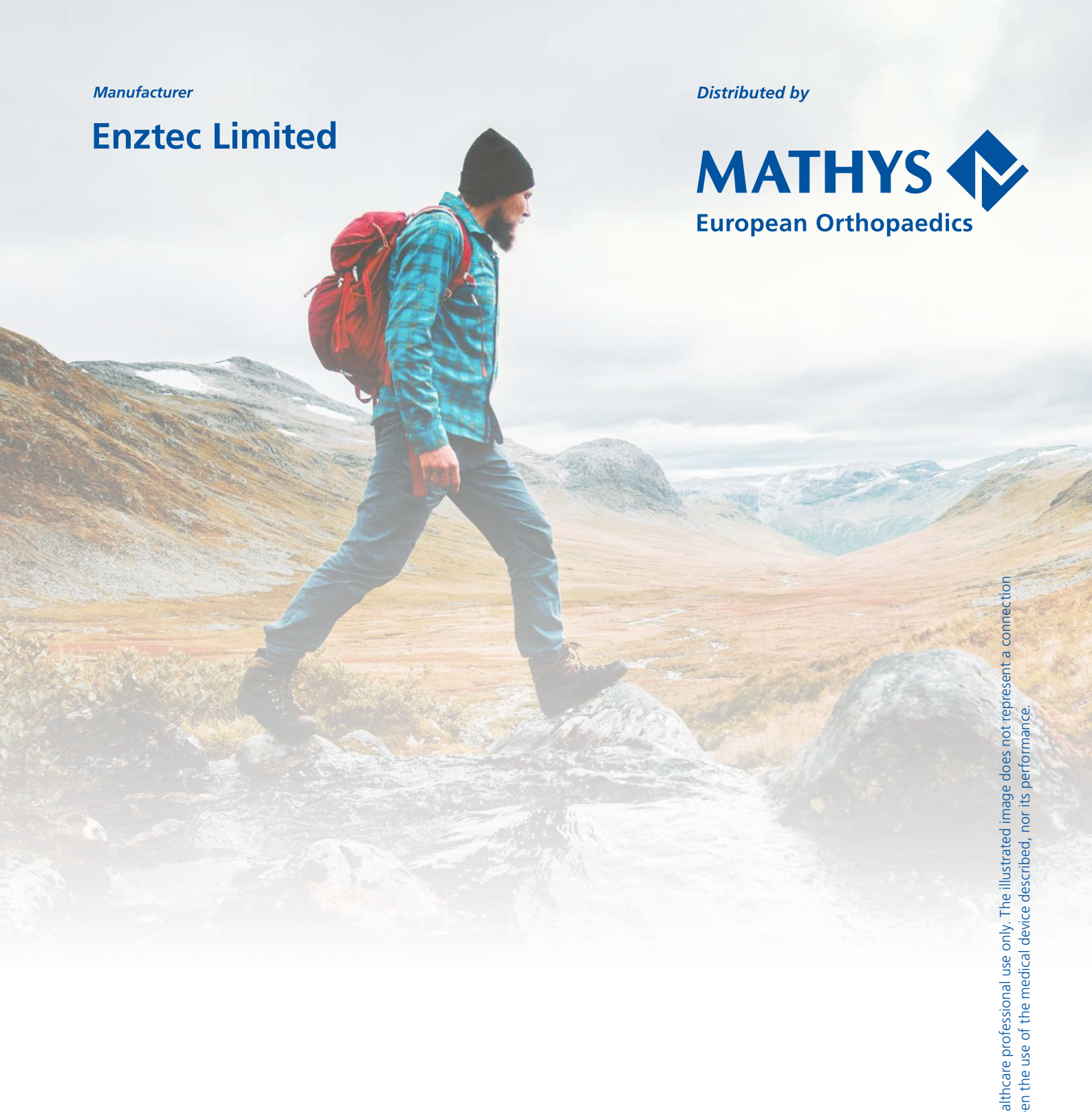

For healthcare professional use only. The illustrated image does not represent a connection between the use of the medical device described, nor its performance.

Preservation in motion



Offset reamer driver

*Building on our heritage
Moving technology forward
Step by step with our clinical partners
Towards a goal of preserving mobility*

Preservation in motion

As a Swiss company, Mathys is committed to this guiding principle and pursues a product portfolio with the goal of further developing traditional philosophies with respect to materials or design in order to address existing clinical challenges. This is reflected in our imagery: traditional Swiss activities in conjunction with continuously evolving sporting equipment.

Offset reamer driver

The offset reamer driver is combined with an acetabular reamer to prepare the acetabulum for cup implantation. The surgeon can **choose from two types of acetabular reamer connection**: locking or open capture. The locking option is the standard version. It also allows reverse reaming if bone chips are to be compressed for defect filling in the acetabulum during bone grafting. The open option allows the user to connect the reamer driver to the acetabular reamer in situ and to disconnect it there as well.


In addition, the **low profile construction of the offset reamer driver** where it connects to the reamer enables good handling regardless of approach. The **handle position** can also be adjusted intraoperatively to help the surgeon achieve the best possible reaming position and force transmission.



Advantages




locking

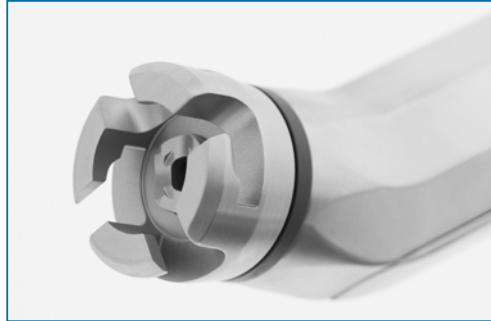

open

Two reamer connections

- Two types of acetabular reamer connection: locking or open capture
- Bone grafting possible by means of reverse reaming with the locking reamer driver
- In situ reamer capture possible with the open reamer driver
- Good handling with all common approaches to the hip due to the low profile design

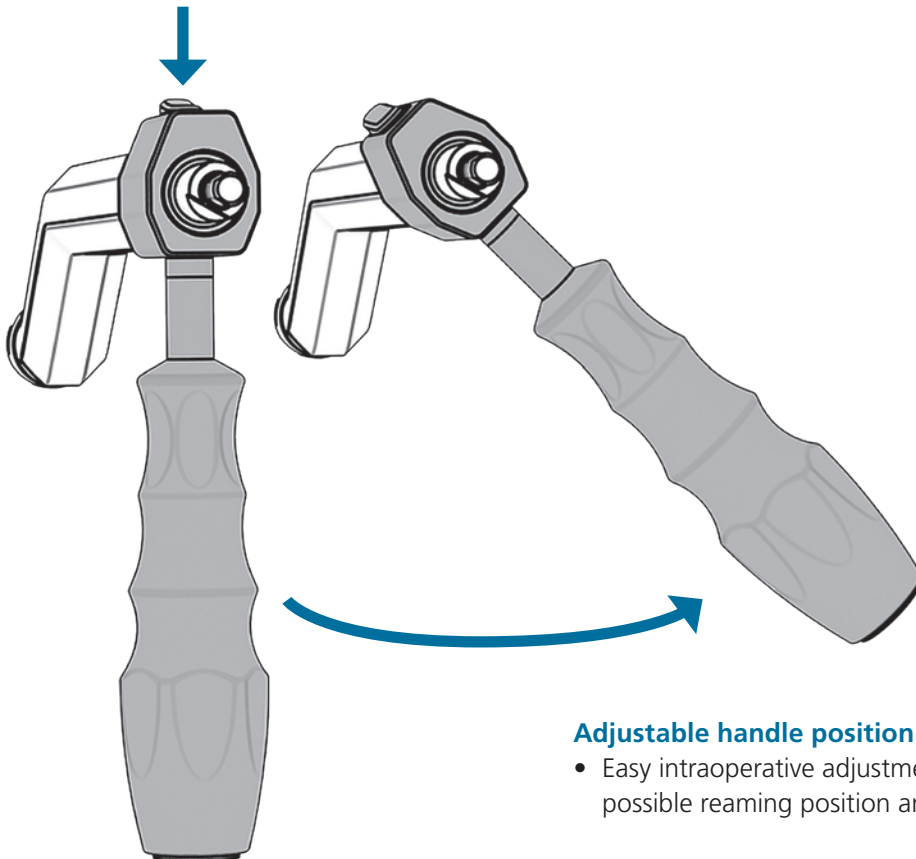

Adjustable handle position

- Easy intraoperative adjustment of the handle allows best possible reaming position and good force transmission.

Modular design

- Intuitive assembly and simple disassembly for cleaning the instrument
- Replacement parts available for repair exchange
- Options available for common power tool connections



Offset body



Offset drive train



Offset coupling



Offset cover



Handle

Instruments



Item no.	Description
51.34.1150A	Offset reamer driver – Locking – AO
51.34.1151A	Offset reamer driver – Open – AO
51.34.1169A	Offset reamer driver – Locking – Zimmer
51.34.1170A	Offset reamer driver – Open – Zimmer
51.34.1171A	Offset reamer driver – Locking – Hudson
51.34.1172A	Offset reamer driver – Open – Hudson

Replacement parts:

Item no.	Description
4250-7048	Offset body
4250-7031	Offset cover
4250-7012	Handle

Item no.	Description
4250-7032	Offset coupling – Zimmer
4250-7033	Offset coupling – Hudson
4250-7034	Offset coupling – Large AO

Item no.	Description
4250-7035	Offset drive train – Open
4250-7036	Offset drive train – Locking – Crossbar

Manufacturer

Enztec Limited

Distributed by



Australia Mathys Orthopaedics Pty Ltd
Artarmon, NSW 2064
Tel: +61 2 9417 9200
info.au@mathysmedical.com

Austria Mathys Orthopädie GmbH
2351 Wiener Neudorf
Tel: +43 2236 860 999
info.at@mathysmedical.com

Belgium Mathys Orthopaedics Belux N.V.-S.A.
3001 Leuven
Tel: +32 16 38 81 20
info.be@mathysmedical.com

France Mathys Orthopédie S.A.S
63360 Gerzat
Tel: +33 4 73 23 95 95
info.fr@mathysmedical.com

Germany Mathys Orthopädie GmbH
«Centre of Excellence Sales» Bochum
44809 Bochum
Tel: +49 234 588 59 0
sales.de@mathysmedical.com

«Centre of Excellence Ceramics» Mörsdorf
07646 Mörsdorf/Thür.
Tel: +49 364 284 94 0
info.de@mathysmedical.com

«Centre of Excellence Production» Hermsdorf
07629 Hermsdorf
Tel: +49 364 284 94 110
info.de@mathysmedical.com

Italy Mathys Ortopedia S.r.l.
20141 Milan
Tel: +39 02 4959 8085
info.it@mathysmedical.com

Japan Mathys KK
Tokyo 108-0075
Tel: +81 3 3474 6900
info.jp@mathysmedical.com

New Zealand Mathys Ltd.
Auckland
Tel: +64 9 478 39 00
info.nz@mathysmedical.com

Netherlands Mathys Orthopaedics B.V.
3001 Leuven
Tel: +31 88 1300 500
info.nl@mathysmedical.com

P. R. China Mathys (Shanghai) Medical Device Trading Co., Ltd
Shanghai, 200041
Tel: +86 21 6170 2655
info.cn@mathysmedical.com

Switzerland Mathys (Schweiz) GmbH
2544 Bettlach
Tel: +41 32 644 1 458
info@mathysmedical.com

United Kingdom Mathys Orthopaedics Ltd
Alton, Hampshire GU34 2QL
Tel: +44 8450 580 938
info.uk@mathysmedical.com

Local Marketing Partners in over 30 countries worldwide ...

Manufactured by Enztec Limited • 3/17 Print Place • Middleton • Christchurch 8024 • New Zealand • www.enztec.com

Distributed by Mathys Ltd Bettlach • Robert Mathys Strasse 5 • P.O. Box • 2544 Bettlach • Switzerland • www.mathysmedical.com