

Portfolio 2023

PRODUCT PORTFOLIO




Enovis Continuum of Care

Enovis[™] is the only company with products across the full spectrum of care, from pain management and improved stability to surgery and rehabilitation. Taken together, the Enovis product range delivers an unparalleled continuum of care.

Our strategy is supported by an internal organization of product-based business units. This allows us to respond quickly to market demands, with the focus on developing a portfolio of cutting-edge products and services.

Our Prevention and Recovery division encompasses non-invasive products and solutions, with a dedicated business unit focused on Bracing and Supports, Recovery Science (medical technology and capital equipment for rehabilitation), and FootCare Solutions (medical shoes and hosiery). Moreover, we have the Surgical division focusing on manufacturing and distributing orthopedic implants and complementary associated technologies.

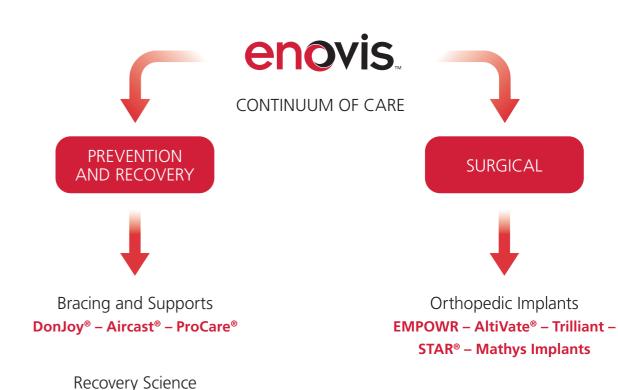

Chattanooga® - CMF - Compex® -

Companion


FootCare Solutions

Dr. Comfort® - Aircast

Instrumentation Adaptable

Navigation Technology

Arvis

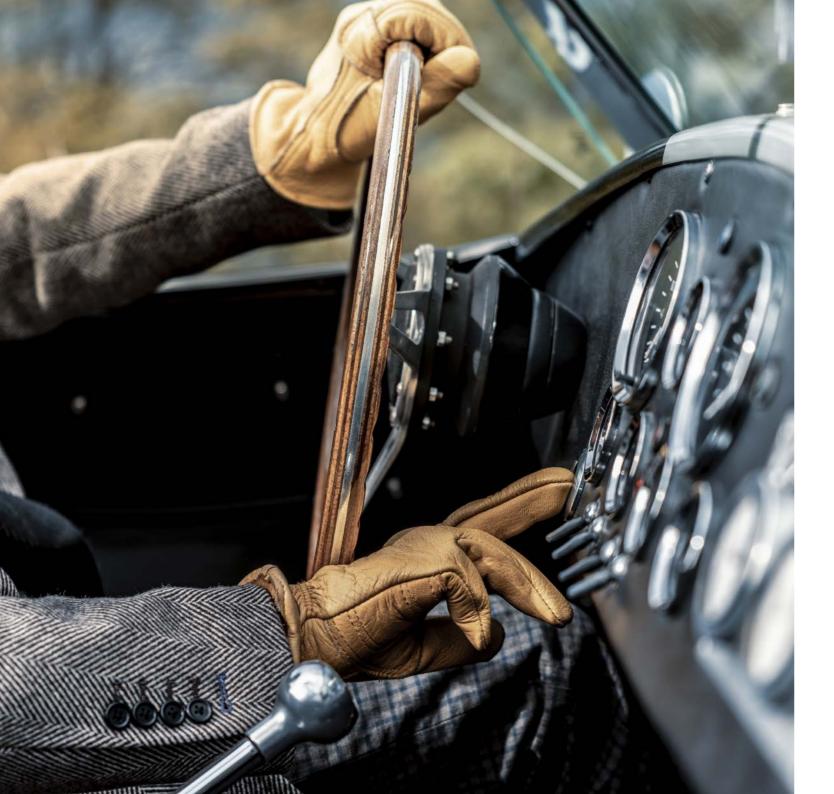

SURGICAL EQUIPMENT

Adaptable® – Hip/Knee/Shoulder





HIP SEGMENT

RM Pressfit vitamys, ceramic head & optimys



Hip primary

Monobloc



RM Pressfit vitamys



RM Pressfit



RM Classic cup

Modular



aneXys Flex/aneXys Uno



seleXys PC



aneXys Cluster/aneXys Multi



seleXys ceramys inlay



aneXys ceramys inlay



aneXys vitamys inlay



seleXys PE inlay

Hip primary

Dual mobility



DS Evolution uncemented





cemented





DS Evolution Revision



DS Evolution PE inlay

Bipolar/hemi



Bipolar



Hemi

Heads



ceramys



symarec



Stainless steel



CoCrMo

Hip primary





Hip revision





KNEE SEGMENT

balanSys BICONDYLAR & EMPOWR Knee System



Partial **knee**

Total knee

balanSys UNI






EMPOWR 3D/PS Knee



Total knee



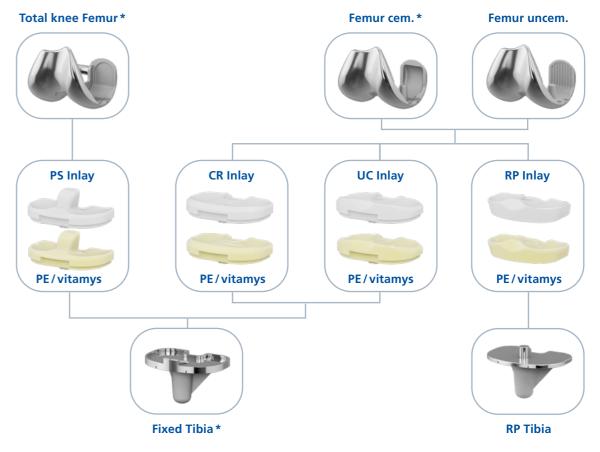

* available in TiNbN













3 Peg

3 Peg FLAT



REVISION knee

balanSys REV





SHOULDER SEGMENT

Affins Short, Affinis Inverse & AltiVate Reverse




Anatomic shoulder

Affinis Short/Affinis Total




Shoulder Fracture / Fracture Reverse

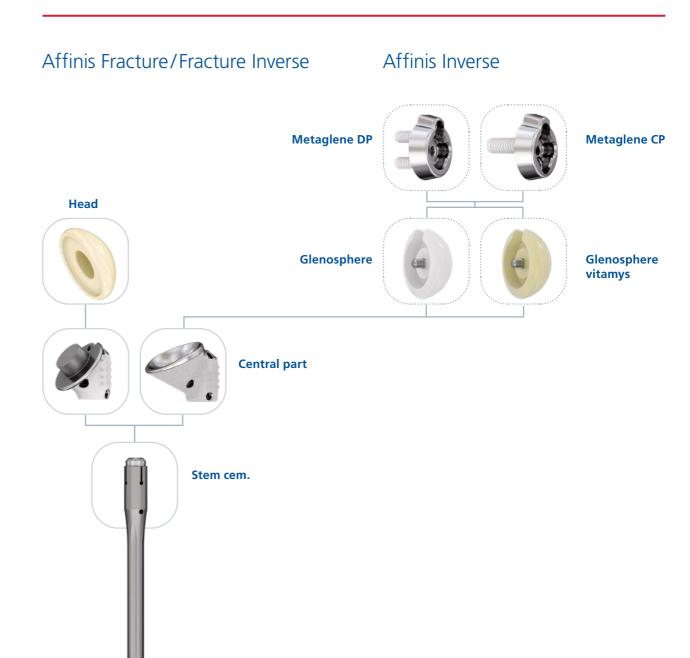



Reverse **shoulder**

Affinis Inverse

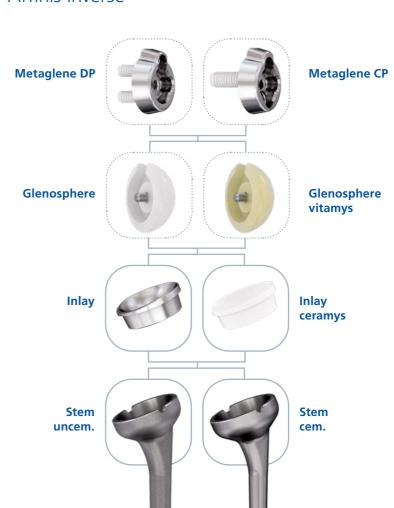

AltiVate® Reverse



Glenoid Baseplate




Glenosphere



Socket Insert

Humeral Stem, standard



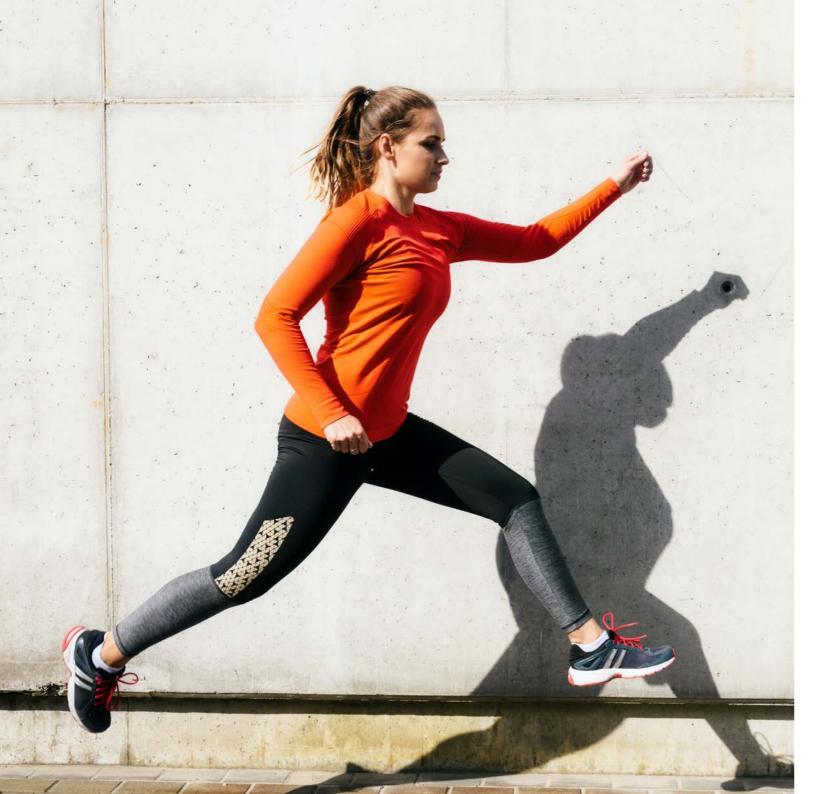

SPORTHOPAEDICS


Ligamys



PROGRESSIVE MATERIALS


vitamys

Mathys was the first manufacturer of orthopaedic prostheses to launch a highly cross-linked polyethylene enriched with vitamin E (VEPE) – a so-called «blended antioxidant highly cross-linked polyethylene (AOHXLPE)» –, trademarked as vitamys. The vitamin-E-enriched, highly crosslinked vitamys is designed to deliver maximum resilience and a long service life. ^{2, 3} Properties such as high resistance to oxidation, ageing and wear, as well as excellent mechanical properties distinguish this proven material and allow long-term survival of the endoprosthesis. ^{2, 3, 4}

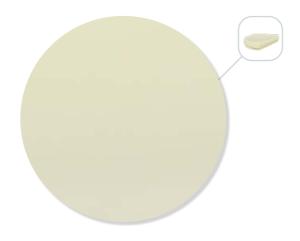


Ceramics

Since the early 1970s, Mathys has been active in the research, development and manufacture of bioceramics, because we are convinced of their advantages: low wear rates, high strength and toughness, good wettability and biologically inert behaviour. ^{1, 5, 6, 7} This makes ceramics a treatment solution not only for young and active patients.

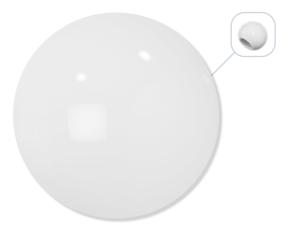

RM titanium coating

Mathys' proprietary titanium coating «RM» (Robert Mathys) enables osseointegration of the implant and permits cementless anchorage. The titanium particles are anchored individually in the polyethylene and not structurally connected to each other. Thus, the elasticity of the implant is not changed by the coating.

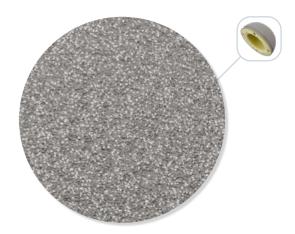



References

- ¹ Lerf R, Wuttke V, Reimelt I, Dallmann F, Delfosse D, Tribological Behaviour of the «Reverse» Inverse Shoulder Prosthesis. 7th International UHMWPE Meeting, Philadelphia 2015.
- Delfosse D, Lerf R, Adlhart C. What happens to the vitamin E in a vitamin-stabilised HXLPE? Karl Knahr (Ed.), Tribology in Total Hip and Knee Arthroplasty. Book Chapter, 2014.
- ³ Lerf R, Zurbrugg D, Delfosse D. Use of vitamin E to protect cross-linked UHMWPE from oxidation. Biomaterials, 2010. 31(13): p. 3643-8.
- ⁴ Bowden AE, Bergström J. Computer Modeling and Simulation of UHMWPE. Kurtz SM (Ed.), UHMWPE Biomaterials Handbook (3rd Edition), 2016.
- ⁵ Willmann G. Improving Bearing Surfaces of Artificial Joints. Adv Eng Mater, 2001. 3(3): p. 135 141.
- ⁶ Barnes DH, Moavenian A, Sharma A, Best SM. Biocompatibility of Ceramics. ASM Handbook (Vol 23), 2012.
- Alexander JJ, Bell SN, Coghlan J, Lerf R, Dallmann F. The effect of vitamin E–enhanced cross-linked polyethylene on wear in shoulder arthroplasty-a wear simulator study. J Shoulder Elbow Surg. 2019 Sep;28(9):1771-1778.



Mathys Ltd Bettlach

Robert Mathys Strasse 5 • 2544 Bettlach • P.O. Box • Switzerland Tel. +41 32 644 1 644 • Fax +41 32 644 1 161 www.mathysmedical.com